

Summary Report of IQC program for G6PD Quantitative Test - Medicon Group - February 2022 -


I. The statistic results of all laboratories in this month

| G6PD | Control N (Lot No.:BJ0922N) | Control D (Lot No.: BJ0922D) |
|-----------------------------|--------------------------------|---------------------------------|
| Labs | 9 | 9 |
| Received results number (n) | 30 | 30 |
| Median | 14.4 (U/g Hb) | 4.7 (U/g Hb) |
| Mean | 14.3 (U/g Hb) | 4.7 (U/g Hb) |
| SD | 0.4 | 0.2 |
| CV | 2.8% | 4.3% |
| Range of G6PD | 13.5 ~ 15.0 (U/g Hb) | 4.3 ~ 4.9 (U/g Hb) |
| Range of Hb | 2.1 ~ 2.6 (g/dL) | 1.9 ~ 2.3 (g/dL) |

*The statistic results are calculated from all labs reported in this month

** G6PD Method = Medicon reagent kit, 37°C

II. The distribution of G6PD reported for each lab in this survey

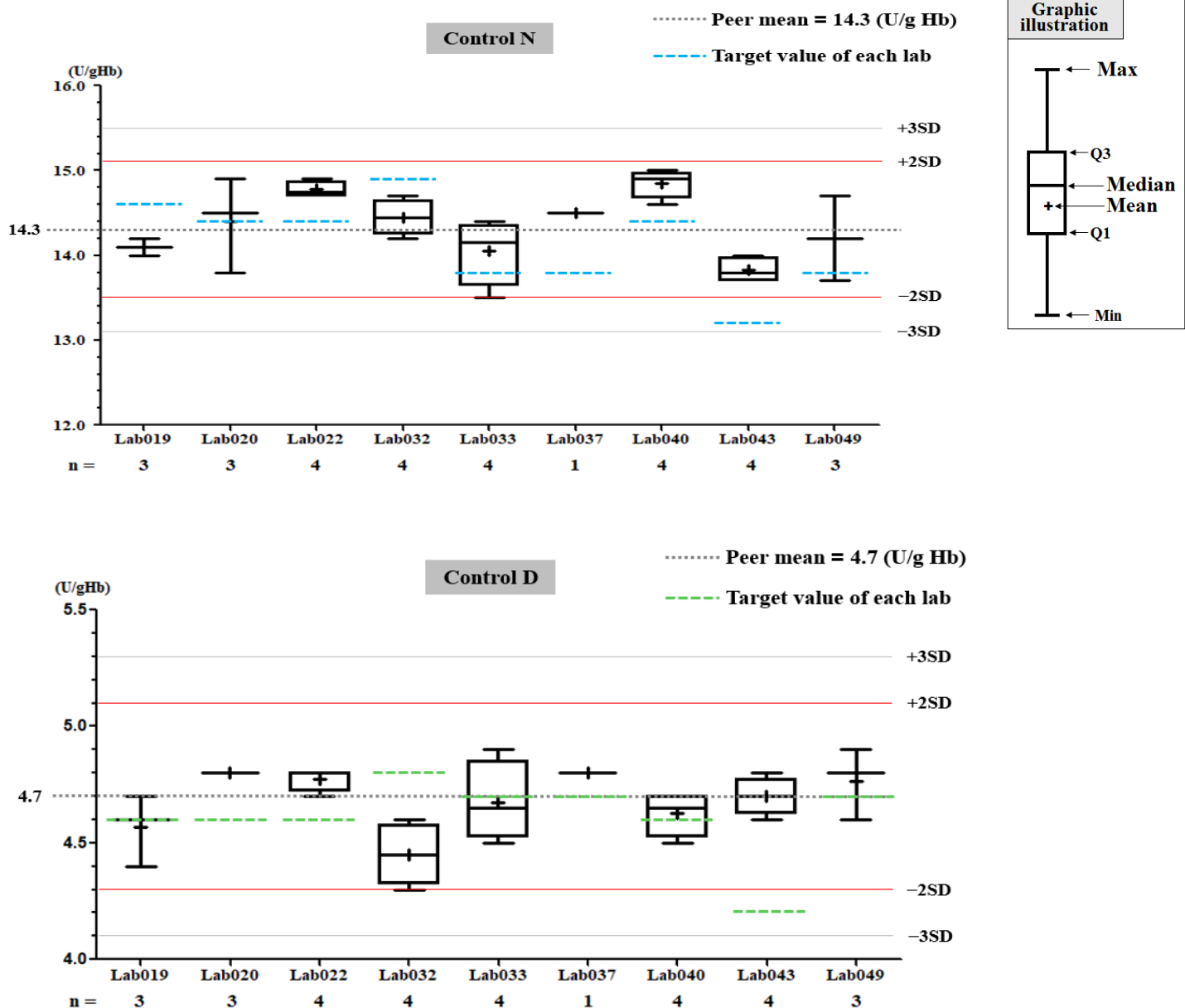

QC Chart of Internal Quality Control (IQC) for G6PD Quantitative Test

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

[Print Table](#)

Lab040

| QC Control Lot No. | Control N | | Control D | |
|---------------------------|-----------------|-----------------------------|-----------------|-----------------------------|
| | BJ0922N | | BJ0922D | |
| Duration of the Analyzing | Month (2023/02) | CUM (2022/09/30~2023/02/28) | Month (2023/02) | CUM (2022/09/30~2023/02/28) |
| Runs (N) | 4 | 26 | 4 | 26 |
| Mean (U/gHb) | 14.9 | 14.8 | 4.6 | 4.7 |
| SD | 0.2 | 0.2 | 0.1 | 0.1 |
| CV (%) | 1.3 | 1.4 | 2.2 | 2.1 |
| Target Value (U/gHb) | 14.4 | 14.4 | 4.6 | 4.6 |
| Total Error (%) | 6.2 | 5.5 | 4.3 | 6.4 |
| TEa (%) | 20 | 20 | 20 | 20 |
| σ | >6 | >6 | >6 | >6 |

Bias (%) = [(| Mean - Target |) / Target] x 100%
 TE : Total Error(%) = Bias (%) + 2 x CV (%)
 σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2023 02 Change ; Cumulative : from 2022 09 30 to 2023 02 28 Change


[\[TOP\]](#)

Control N and Control D SDI QC Chart

Lot No.: BJ0922N ; Duration : 2022-09-30 to 2023-02-28 ; Target : 14.4 ; SD : 0.80 (●)

Lot No.: BJ0922D ; Duration : 2022-09-30 to 2023-02-28 ; Target : 4.6 ; SD : 0.20 (□)

Lab040

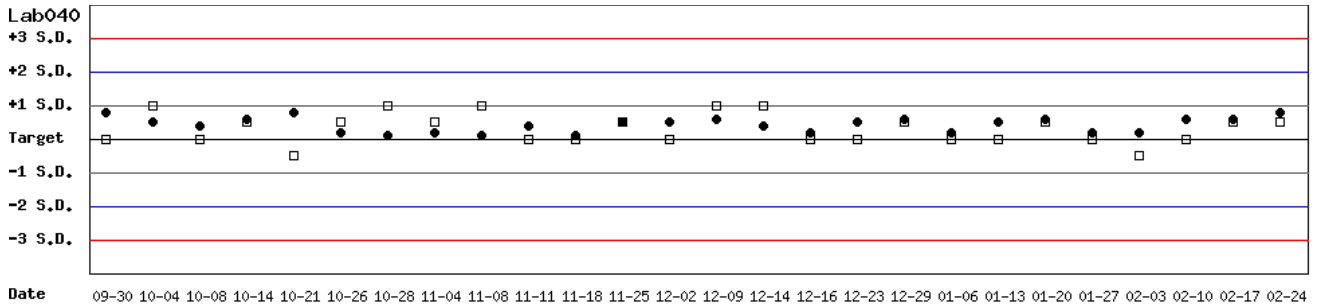


Month : 2023 02 Change ; Cumulative : from 2022 09 30 to 2023 02 28 Change

[\[TOP\]](#)

Control N SDI QC Chart

Lot No.: BJ0922N ; Duration : 2022-09-30 to 2023-02-28 ; Target : 14.4 ; SD : 0.80

Lab040

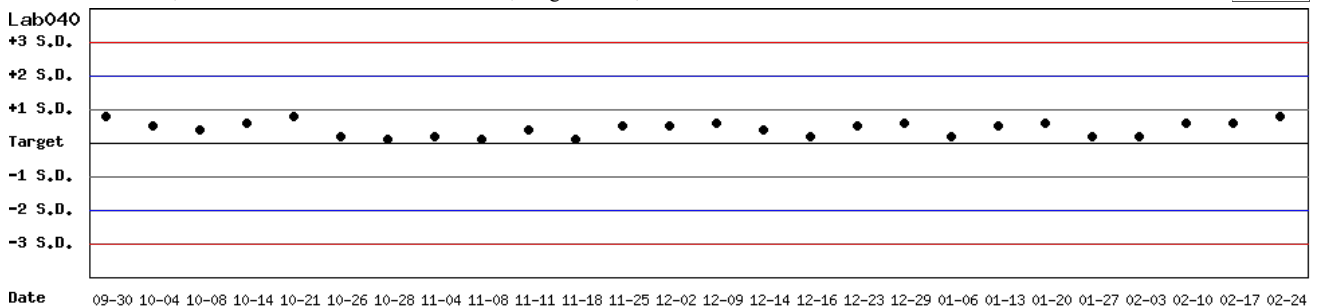


Month : 2023 02 Change ; Cumulative : from 2022 09 30 to 2023 02 28 Change

[\[TOP\]](#)

Control D SDI QC Chart

Lab040

Lot No.: BJ0922D ; Duration : 2022-09-30 to 2023-02-28 ; Target : 4.6 ; SD : 0.20

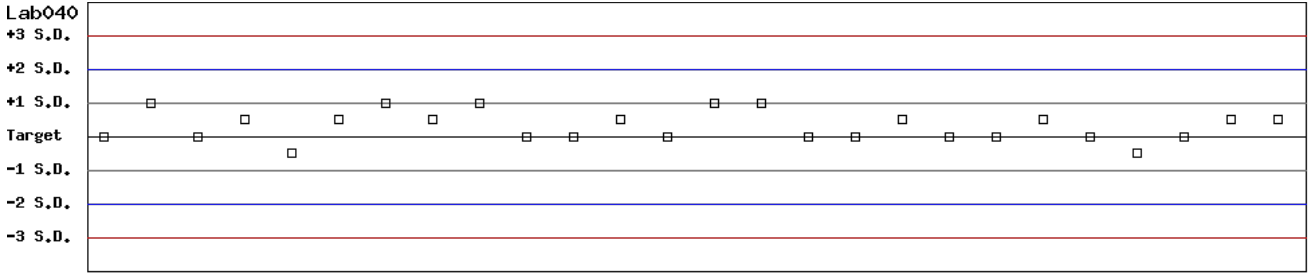

Date 09-30 10-04 10-08 10-14 10-21 10-26 10-28 11-04 11-08 11-11 11-18 11-25 12-02 12-09 12-14 12-16 12-23 12-29 01-06 01-13 01-20 01-27 02-03 02-10 02-17 02-24

Month : ; Cumulative : from to

[TOP](#)

Peer Group Statistics (Table 1)

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

[Print Table 1](#)

Monthly

Month : 2023 02 Change

| UnitID ↑ | Reagent Kit (Code) ↑ | Control N (Lot No.: BJ0922N) | | | | | | | | Control D (Lot No.: BJ0922D) | | | | | | | |
|------------------------|-----------------------------------|------------------------------|--------------|------------|-----|--------|--------|---------|-----|------------------------------|--------------|------------|-----|--------|--------|---------|-----|
| | | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab019 | 5 | 14.6 | 14.1 | 3 | 0.1 | 0.7 | 4.8 | 20 | >6 | 4.6 | 4.6 | 3 | 0.2 | 4.3 | 8.7 | 20 | 4.7 |
| Lab020 | 5 | 14.4 | 14.4 | 3 | 0.6 | 4.2 | 8.3 | 20 | 4.8 | 4.6 | 4.8 | 3 | 0.0 | 0.0 | 4.3 | 20 | >6 |
| Lab022 | 5 | 14.4 | 14.8 | 4 | 0.1 | 0.7 | 4.1 | 20 | >6 | 4.6 | 4.8 | 4 | 0.1 | 2.1 | 8.5 | 20 | >6 |
| Lab032 | 5 | 14.9 | 14.5 | 4 | 0.2 | 1.4 | 5.4 | 20 | >6 | 4.8 | 4.5 | 4 | 0.1 | 2.2 | 10.7 | 20 | >6 |
| Lab033 | 5 | 13.8 | 14.1 | 4 | 0.4 | 2.8 | 7.8 | 20 | >6 | 4.7 | 4.7 | 4 | 0.2 | 4.3 | 8.5 | 20 | 4.7 |
| Lab037 | 5 | 13.8 | 14.5 | 1 | - | - | - | 20 | - | 4.7 | 4.8 | 1 | - | - | - | 20 | - |
| Lab040 | 5 | 14.4 | 14.9 | 4 | 0.2 | 1.3 | 6.2 | 20 | >6 | 4.6 | 4.6 | 4 | 0.1 | 2.2 | 4.3 | 20 | >6 |
| Lab043 | 5 | 13.2 | 13.8 | 4 | 0.2 | 1.4 | 7.4 | 20 | >6 | 4.2 | 4.7 | 4 | 0.1 | 2.1 | 16.2 | 20 | 3.9 |
| Lab049 | 5 | 13.8 | 14.2 | 3 | 0.5 | 3.5 | 9.9 | 20 | 4.9 | 4.7 | 4.8 | 3 | 0.2 | 4.2 | 10.5 | 20 | 4.3 |
| Total | - | - | 14.3 | 30 | 0.4 | 2.8 | - | - | - | - | 4.7 | 30 | 0.2 | 4.3 | - | - | - |

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[TOP](#)

Cumulative

Cumulative : from 2021 01 01 to 2023 02 28 Change

| UnitID ↑ | Reagent Kit (Code) ↑ | Control N (Lot No.: BJ0922N) | | | | | | | | Control D (Lot No.: BJ0922D) | | | | | | | |
|------------------------|-----------------------------------|------------------------------|--------------|------------|-----|--------|--------|---------|-----|------------------------------|--------------|------------|-----|--------|--------|---------|-----|
| | | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab018 | 5 | 16.2 | 16.4 | 8 | 0.4 | 2.4 | 6.1 | 20 | >6 | 5.5 | 5.5 | 8 | 0.2 | 3.6 | 7.3 | 20 | 5.6 |
| Lab019 | 5 | 14.6 | 14.3 | 47 | 0.5 | 3.5 | 9.0 | 20 | 5.1 | 4.6 | 4.6 | 47 | 0.2 | 4.3 | 8.7 | 20 | 4.7 |
| Lab020 | 5 | 14.4 | 14.4 | 55 | 0.5 | 3.5 | 6.9 | 20 | 5.7 | 4.6 | 4.7 | 55 | 0.1 | 2.1 | 6.4 | 20 | >6 |
| Lab022 | 5 | 14.4 | 14.4 | 21 | 0.5 | 3.5 | 6.9 | 20 | 5.7 | 4.6 | 4.6 | 21 | 0.2 | 4.3 | 8.7 | 20 | 4.7 |
| Lab028 | 5 | 15.2 | 15.1 | 34 | 1.1 | 7.3 | 15.2 | 20 | 2.6 | 4.9 | 4.9 | 34 | 0.3 | 6.1 | 12.2 | 20 | 3.3 |
| Lab032 | 5 | 14.9 | 14.4 | 114 | 0.4 | 2.8 | 8.9 | 20 | 5.9 | 4.8 | 4.5 | 114 | 0.1 | 2.2 | 10.7 | 20 | >6 |
| Lab033 | 5 | 13.8 | 13.7 | 59 | 0.6 | 4.4 | 9.5 | 20 | 4.4 | 4.7 | 4.6 | 59 | 0.3 | 6.5 | 15.2 | 20 | 2.7 |
| Lab037 | 5 | 13.8 | 14.6 | 49 | 0.4 | 2.7 | 11.3 | 20 | 5.3 | 4.7 | 4.7 | 49 | 0.2 | 4.3 | 8.5 | 20 | 4.7 |
| Lab040 | 5 | 14.4 | 14.7 | 39 | 0.2 | 1.4 | 4.8 | 20 | >6 | 4.6 | 4.6 | 39 | 0.1 | 2.2 | 4.3 | 20 | >6 |
| Lab043 | 5 | 13.2 | 13.7 | 77 | 1.4 | 10.2 | 24.2 | 20 | 1.6 | 4.2 | 4.4 | 77 | 0.5 | 11.4 | 27.5 | 20 | 1.3 |
| Lab049 | 5 | 13.8 | 14.6 | 87 | 0.7 | 4.8 | 15.4 | 20 | 3.0 | 4.7 | 4.7 | 87 | 0.2 | 4.3 | 8.5 | 20 | 4.7 |
| Total | - | - | 14.3 | 590 | 0.9 | 6.3 | - | - | - | - | 4.6 | 590 | 0.3 | 6.5 | - | - | - |

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[TOP](#)

| Reagent Kit | Reagent Code |
|-------------|--------------|
| Medicon | 5 |

Peer Group Statistics (Table 2)

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

[Print Table 2](#)

Control N Month vs. Cumulative

| | | Control N (Lot No.: BJ0922N) | | | | | | | | | | | | | | | |
|------------------------|----------------------|------------------------------|--------------|------------|-----|--------|--------|---------|-----|-----------------------------|--------------|------------|-----|--------|--------|---------|-----|
| | | Month (2023/02) | | | | | | | | CUM (2021/01/01~2023/02/28) | | | | | | | |
| UnitID ↑ | Reagent Kit (Code) ↑ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab018 | 5 | 16.2 | - | 0 | - | - | - | 20 | - | 16.2 | 16.4 | 8 | 0.4 | 2.4 | 6.1 | 20 | >6 |
| Lab019 | 5 | 14.6 | 14.1 | 3 | 0.1 | 0.7 | 4.8 | 20 | >6 | 14.6 | 14.3 | 47 | 0.5 | 3.5 | 9.0 | 20 | 5.1 |
| Lab020 | 5 | 14.4 | 14.4 | 3 | 0.6 | 4.2 | 8.3 | 20 | 4.8 | 14.4 | 14.4 | 55 | 0.5 | 3.5 | 6.9 | 20 | 5.7 |
| Lab022 | 5 | 14.4 | 14.8 | 4 | 0.1 | 0.7 | 4.1 | 20 | >6 | 14.4 | 14.4 | 21 | 0.5 | 3.5 | 6.9 | 20 | 5.7 |
| Lab028 | 5 | 15.2 | - | 0 | - | - | - | 20 | - | 15.2 | 15.1 | 34 | 1.1 | 7.3 | 15.2 | 20 | 2.6 |
| Lab032 | 5 | 14.9 | 14.5 | 4 | 0.2 | 1.4 | 5.4 | 20 | >6 | 14.9 | 14.4 | 114 | 0.4 | 2.8 | 8.9 | 20 | 5.9 |
| Lab033 | 5 | 13.8 | 14.1 | 4 | 0.4 | 2.8 | 7.8 | 20 | >6 | 13.8 | 13.7 | 59 | 0.6 | 4.4 | 9.5 | 20 | 4.4 |
| Lab037 | 5 | 13.8 | 14.5 | 1 | - | - | - | 20 | - | 13.8 | 14.6 | 49 | 0.4 | 2.7 | 11.3 | 20 | 5.3 |
| Lab040 | 5 | 14.4 | 14.9 | 4 | 0.2 | 1.3 | 6.2 | 20 | >6 | 14.4 | 14.7 | 39 | 0.2 | 1.4 | 4.8 | 20 | >6 |
| Lab043 | 5 | 13.2 | 13.8 | 4 | 0.2 | 1.4 | 7.4 | 20 | >6 | 13.2 | 13.7 | 77 | 1.4 | 10.2 | 24.2 | 20 | 1.6 |
| Lab049 | 5 | 13.8 | 14.2 | 3 | 0.5 | 3.5 | 9.9 | 20 | 4.9 | 13.8 | 14.6 | 87 | 0.7 | 4.8 | 15.4 | 20 | 3.0 |
| Total | - | - | 14.3 | 30 | 0.4 | 2.8 | - | - | - | - | 14.3 | 590 | 0.9 | 6.3 | - | - | - |

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2023 02 Change

Cumulative : from 2021 01 01 to 2023 02 28 Change

[TOP](#)

Control D Month vs. Cumulative

| | | Control D (Lot No.: BJ0922D) | | | | | | | | | | | | | | | |
|------------------------|----------------------|------------------------------|--------------|------------|-----|--------|--------|---------|-----|-----------------------------|--------------|------------|-----|--------|--------|---------|-----|
| | | Month (2023/02) | | | | | | | | CUM (2021/01/01~2023/02/28) | | | | | | | |
| UnitID ↑ | Reagent Kit (Code) ↑ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab018 | 5 | 5.5 | - | 0 | - | - | - | 20 | - | 5.5 | 5.5 | 8 | 0.2 | 3.6 | 7.3 | 20 | 5.6 |
| Lab019 | 5 | 4.6 | 4.6 | 3 | 0.2 | 4.3 | 8.7 | 20 | 4.7 | 4.6 | 4.6 | 47 | 0.2 | 4.3 | 8.7 | 20 | 4.7 |
| Lab020 | 5 | 4.6 | 4.8 | 3 | 0.0 | 0.0 | 4.3 | 20 | >6 | 4.6 | 4.7 | 55 | 0.1 | 2.1 | 6.4 | 20 | >6 |
| Lab022 | 5 | 4.6 | 4.8 | 4 | 0.1 | 2.1 | 8.5 | 20 | >6 | 4.6 | 4.6 | 21 | 0.2 | 4.3 | 8.7 | 20 | 4.7 |
| Lab028 | 5 | 4.9 | - | 0 | - | - | - | 20 | - | 4.9 | 4.9 | 34 | 0.3 | 6.1 | 12.2 | 20 | 3.3 |
| Lab032 | 5 | 4.8 | 4.5 | 4 | 0.1 | 2.2 | 10.7 | 20 | >6 | 4.8 | 4.5 | 114 | 0.1 | 2.2 | 10.7 | 20 | >6 |
| Lab033 | 5 | 4.7 | 4.7 | 4 | 0.2 | 4.3 | 8.5 | 20 | 4.7 | 4.7 | 4.6 | 59 | 0.3 | 6.5 | 15.2 | 20 | 2.7 |
| Lab037 | 5 | 4.7 | 4.8 | 1 | - | - | - | 20 | - | 4.7 | 4.7 | 49 | 0.2 | 4.3 | 8.5 | 20 | 4.7 |
| Lab040 | 5 | 4.6 | 4.6 | 4 | 0.1 | 2.2 | 4.3 | 20 | >6 | 4.6 | 4.6 | 39 | 0.1 | 2.2 | 4.3 | 20 | >6 |
| Lab043 | 5 | 4.2 | 4.7 | 4 | 0.1 | 2.1 | 16.2 | 20 | 3.9 | 4.2 | 4.4 | 77 | 0.5 | 11.4 | 27.5 | 20 | 1.3 |
| Lab049 | 5 | 4.7 | 4.8 | 3 | 0.2 | 4.2 | 10.5 | 20 | 4.3 | 4.7 | 4.7 | 87 | 0.2 | 4.3 | 8.5 | 20 | 4.7 |
| Total | - | - | 4.7 | 30 | 0.2 | 4.3 | - | - | - | - | 4.6 | 590 | 0.3 | 6.5 | - | - | - |

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2023 02 Change

Cumulative : from 2021 01 01 to 2023 02 28 Change

[TOP](#)

| Reagent Kit | Reagent Code |
|-------------|--------------|
| Medicon | 5 |