

Summary Report of IQC program for G6PD Quantitative Test - Medicon Group - April 2023 -


I. The statistic results of all laboratories in this month

G6PD	Control N (Lot No.:BJ0922N)	Control D (Lot No.: BJ0922D)
Labs	9	9
Received results number (n)	30	30
Median	14.3 (U/g Hb)	4.7 (U/g Hb)
Mean	14.1 (U/g Hb)	4.7 (U/g Hb)
SD	0.5	0.3
CV	3.5%	4.3%
Range of G6PD	12.8 ~ 14.9 (U/g Hb)	4.3 ~ 5.3 (U/g Hb)
Range of Hb	2.1 ~ 2.8 (g/dL)	1.9 ~ 2.4 (g/dL)

*The statistic results are calculated from all labs reported in this month

** G6PD Method = Medicon reagent kit, 37°C

II. The distribution of G6PD reported for each lab in this survey

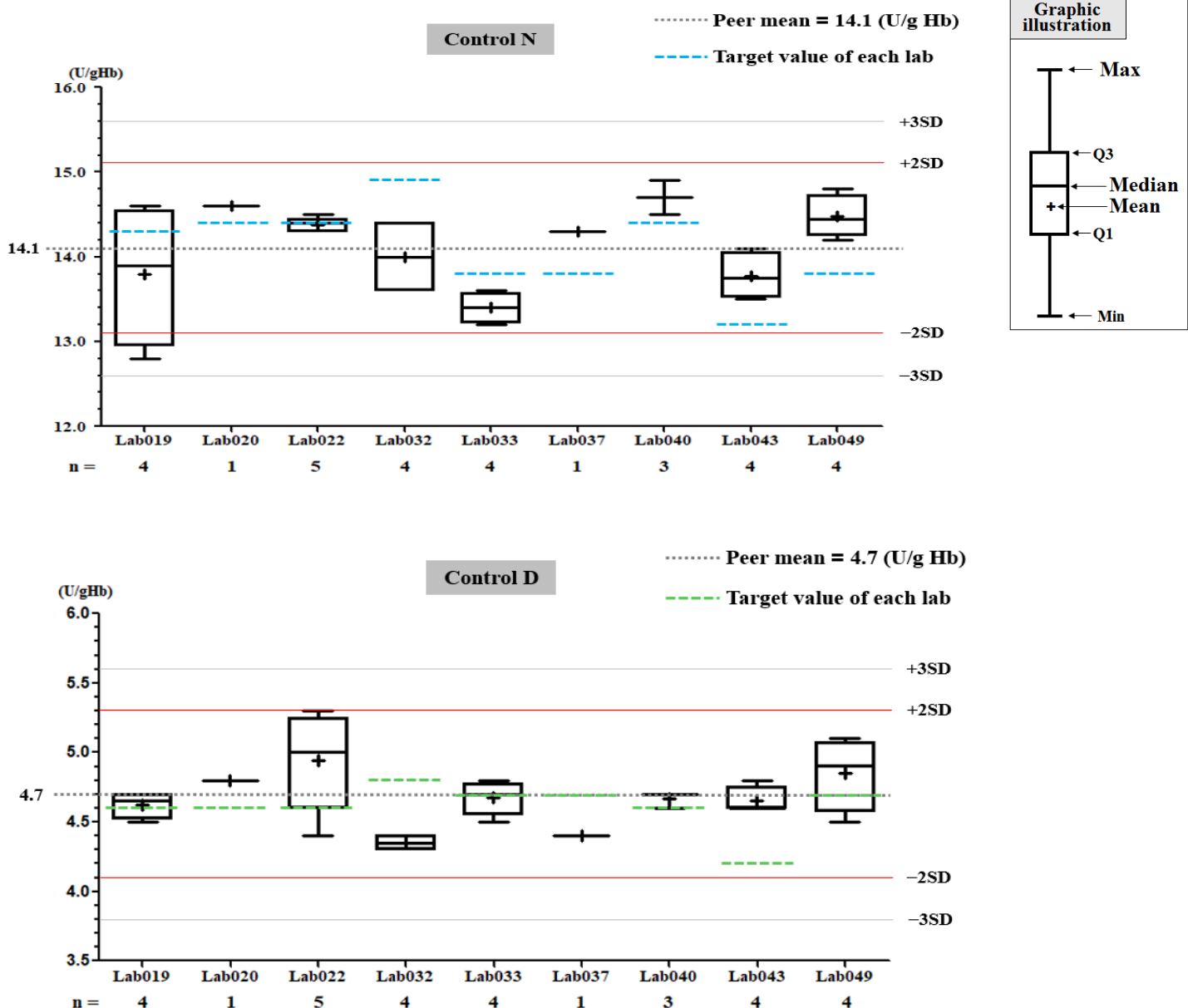

QC Chart of Internal Quality Control (IQC) for G6PD Quantitative Test

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Print Table

Lab037

QC Control Lot No.	Control N		Control D	
	BJ0922N		BJ0922D	
Duration of the Analyzing	Month (2023/04)	CUM (2022/06/16~2023/04/30)	Month (2023/04)	CUM (2022/06/16~2023/04/30)
Runs (N)	1	26	1	26
Mean (U/gHb)	14.3	14.3	4.4	4.7
SD	-	0.4	-	0.2
CV (%)	-	2.8	-	4.3
Target Value (U/gHb)	13.8	13.8	4.7	4.7
Total Error (%)	-	9.2	-	8.5
TEa (%)	20	20	20	20
σ	-	5.8	-	4.7

Bias (%) = [(| Mean - Target |) / Target] x 100%
 TE : Total Error(%) = Bias (%) + 2 x CV (%)
 σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2023 04 Change ; Cumulative : from 2022 06 16 to 2023 04 30 Change


[\[TOP\]](#)

Control N and Control D SDI QC Chart

Lot No.: BJ0922N ; Duration : 2022-06-16 to 2023-04-30 ; Target : 13.8 ; SD : 1.38 (●)

Lot No.: BJ0922D ; Duration : 2022-06-16 to 2023-04-30 ; Target : 4.7 ; SD : 0.47 (□)

Lab037

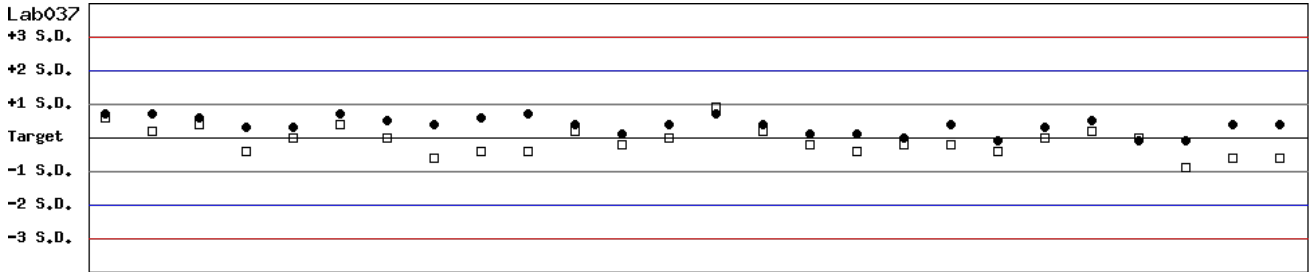

Date 06-16 07-01 07-07 07-14 07-22 07-28 08-11 08-18 08-25 08-31 08-31 09-16 09-29 10-13 10-27 11-10 11-24 12-08 12-22 01-05 01-19 02-02 03-02 03-17 03-30 04-27


Month : 2023 04 Change ; Cumulative : from 2022 06 16 to 2023 04 30 Change

[\[TOP\]](#)

Control N SDI QC Chart

Lot No.: BJ0922N ; Duration : 2022-06-16 to 2023-04-30 ; Target : 13.8 ; SD : 1.38

Lab037

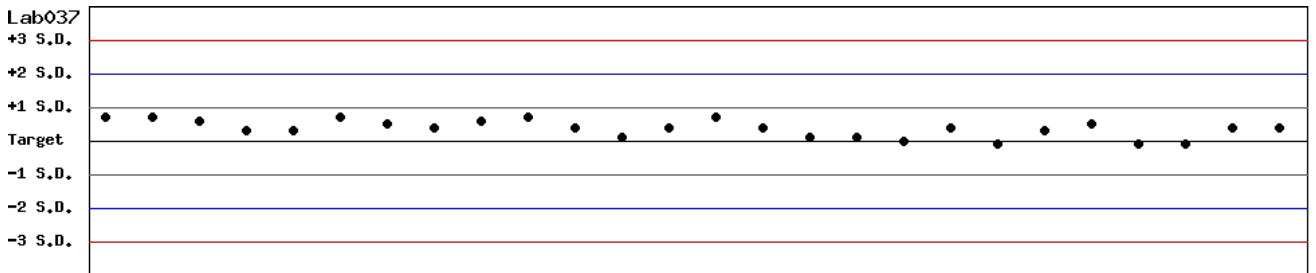

Date 06-16 07-01 07-07 07-14 07-22 07-28 08-11 08-18 08-25 08-31 08-31 09-16 09-29 10-13 10-27 11-10 11-24 12-08 12-22 01-05 01-19 02-02 03-02 03-17 03-30 04-27


Month : 2023 04 Change ; Cumulative : from 2022 06 16 to 2023 04 30 Change

[\[TOP\]](#)

Control D SDI QC Chart

Lab037

Lot No.: BJ0922D ; Duration : 2022-06-16 to 2023-04-30 ; Target : 4.7 ; SD : 0.47

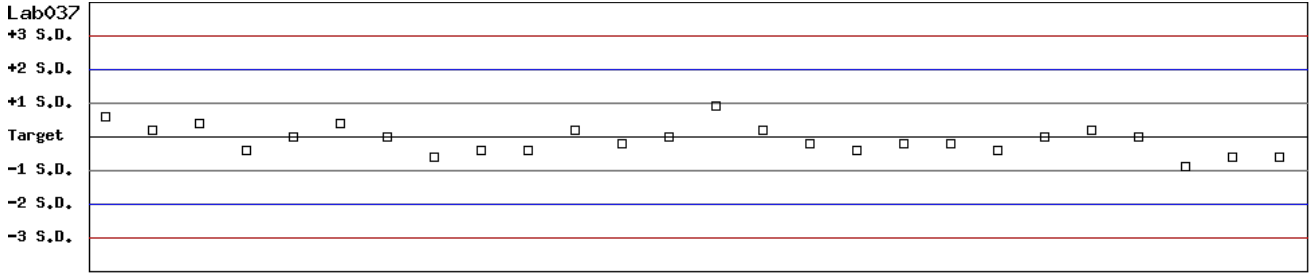

Date 06-16 07-01 07-07 07-14 07-22 07-28 08-11 08-18 08-25 08-31 08-31 09-16 09-29 10-13 10-27 11-10 11-24 12-08 12-22 01-05 01-19 02-02 03-02 03-17 03-30 04-27

Month : ; Cumulative : from to

[TOP](#)

Peer Group Statistics (Table 1)

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

[Print Table 1](#)

Monthly

Month : 2023 04 Change

UnitID ↑	Reagent Kit (Code) ↑	Control N (Lot No.: BJ0922N)								Control D (Lot No.: BJ0922D)							
		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab019	5	14.3	13.8	4	0.8	5.8	15.1	20	2.8	4.6	4.6	4	0.1	2.2	4.3	20	>6
Lab020	5	14.4	14.6	1	-	-	-	20	-	4.6	4.8	1	-	-	-	20	-
Lab022	5	14.4	14.4	5	0.1	0.7	1.4	20	>6	4.6	4.9	5	0.4	8.2	22.8	20	1.6
Lab032	5	14.9	14.0	4	0.5	3.6	13.2	20	3.9	4.8	4.4	4	0.1	2.3	12.9	20	5.1
Lab033	5	13.8	13.4	4	0.2	1.5	5.9	20	>6	4.7	4.7	4	0.1	2.1	4.3	20	>6
Lab037	5	13.8	14.3	1	-	-	-	20	-	4.7	4.4	1	-	-	-	20	-
Lab040	5	14.4	14.7	3	0.2	1.4	4.8	20	>6	4.6	4.7	3	0.1	2.1	6.4	20	>6
Lab043	5	13.2	13.8	4	0.3	2.2	8.9	20	>6	4.2	4.7	4	0.1	2.1	16.2	20	3.9
Lab049	5	13.8	14.5	4	0.3	2.1	9.2	20	>6	4.7	4.9	4	0.3	6.1	16.5	20	2.6
Total	-	-	14.1	30	0.5	3.5	-	-	-	-	4.7	30	0.3	6.4	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Cumulative

Cumulative : from 2021 01 01 to 2023 04 30 Change

UnitID ↑	Reagent Kit (Code) ↑	Control N (Lot No.: BJ0922N)								Control D (Lot No.: BJ0922D)							
		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	16.2	16.4	8	0.4	2.4	6.1	20	>6	5.5	5.5	8	0.2	3.6	7.3	20	5.6
Lab019	5	14.3	14.3	54	0.6	4.2	8.4	20	4.8	4.6	4.6	54	0.2	4.3	8.7	20	4.7
Lab020	5	14.4	14.4	60	0.5	3.5	6.9	20	5.7	4.6	4.7	60	0.1	2.1	6.4	20	>6
Lab022	5	14.4	14.4	27	0.5	3.5	6.9	20	5.7	4.6	4.7	27	0.3	6.4	14.9	20	2.8
Lab028	5	15.2	15.1	34	1.1	7.3	15.2	20	2.6	4.9	4.9	34	0.3	6.1	12.2	20	3.3
Lab032	5	14.9	14.3	123	0.4	2.8	9.6	20	5.7	4.8	4.5	123	0.2	4.4	15.1	20	3.1
Lab033	5	13.8	13.7	68	0.6	4.4	9.5	20	4.4	4.7	4.6	68	0.3	6.5	15.2	20	2.7
Lab037	5	13.8	14.5	53	0.5	3.4	12.0	20	4.4	4.7	4.7	53	0.2	4.3	8.5	20	4.7
Lab040	5	14.4	14.7	47	0.2	1.4	4.8	20	>6	4.6	4.6	47	0.1	2.2	4.3	20	>6
Lab043	5	13.2	13.7	85	1.3	9.5	22.8	20	1.7	4.2	4.5	85	0.5	11.1	29.4	20	1.2
Lab049	5	13.8	14.5	91	0.7	4.8	14.7	20	3.1	4.7	4.7	91	0.2	4.3	8.5	20	4.7
Total	-	-	14.3	650	0.8	5.6	-	-	-	-	4.6	650	0.3	6.5	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Medicon	5

Peer Group Statistics (Table 2)

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) ▼ Change

Select Reagent Kit : 5 - Medicon ▼ Change

[Print Table 2](#)

Control N Month vs. Cumulative

		Control N (Lot No.: BJ0922N)															
		Month (2023/04)								CUM (2021/01/01~2023/04/30)							
UnitID ↑	Reagent Kit (Code) ↑	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	16.2	-	0	-	-	-	20	-	16.2	16.4	8	0.4	2.4	6.1	20	>6
Lab019	5	14.3	13.8	4	0.8	5.8	15.1	20	2.8	14.3	14.3	54	0.6	4.2	8.4	20	4.8
Lab020	5	14.4	14.6	1	-	-	-	20	-	14.4	14.4	60	0.5	3.5	6.9	20	5.7
Lab022	5	14.4	14.4	5	0.1	0.7	1.4	20	>6	14.4	14.4	27	0.5	3.5	6.9	20	5.7
Lab028	5	15.2	-	0	-	-	-	20	-	15.2	15.1	34	1.1	7.3	15.2	20	2.6
Lab032	5	14.9	14.0	4	0.5	3.6	13.2	20	3.9	14.9	14.3	123	0.4	2.8	9.6	20	5.7
Lab033	5	13.8	13.4	4	0.2	1.5	5.9	20	>6	13.8	13.7	68	0.6	4.4	9.5	20	4.4
Lab037	5	13.8	14.3	1	-	-	-	20	-	13.8	14.5	53	0.5	3.4	12.0	20	4.4
Lab040	5	14.4	14.7	3	0.2	1.4	4.8	20	>6	14.4	14.7	47	0.2	1.4	4.8	20	>6
Lab043	5	13.2	13.8	4	0.3	2.2	8.9	20	>6	13.2	13.7	85	1.3	9.5	22.8	20	1.7
Lab049	5	13.8	14.5	4	0.3	2.1	9.2	20	>6	13.8	14.5	91	0.7	4.8	14.7	20	3.1
Total	-	-	14.1	30	0.5	3.5	-	-	-	-	14.3	650	0.8	5.6	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2023 ▼ 04 ▼ Change

Cumulative : from 2021 ▼ 01 ▼ 01 ▼ to 2023 ▼ 04 ▼ 30 ▼ Change

[TOP](#)

Control D Month vs. Cumulative

		Control D (Lot No.: BJ0922D)															
		Month (2023/04)								CUM (2021/01/01~2023/04/30)							
UnitID ↑	Reagent Kit (Code) ↑	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	5.5	-	0	-	-	-	20	-	5.5	5.5	8	0.2	3.6	7.3	20	5.6
Lab019	5	4.6	4.6	4	0.1	2.2	4.3	20	>6	4.6	4.6	54	0.2	4.3	8.7	20	4.7
Lab020	5	4.6	4.8	1	-	-	-	20	-	4.6	4.7	60	0.1	2.1	6.4	20	>6
Lab022	5	4.6	4.9	5	0.4	8.2	22.8	20	1.6	4.6	4.7	27	0.3	6.4	14.9	20	2.8
Lab028	5	4.9	-	0	-	-	-	20	-	4.9	4.9	34	0.3	6.1	12.2	20	3.3
Lab032	5	4.8	4.4	4	0.1	2.3	12.9	20	5.1	4.8	4.5	123	0.2	4.4	15.1	20	3.1
Lab033	5	4.7	4.7	4	0.1	2.1	4.3	20	>6	4.7	4.6	68	0.3	6.5	15.2	20	2.7
Lab037	5	4.7	4.4	1	-	-	-	20	-	4.7	4.7	53	0.2	4.3	8.5	20	4.7
Lab040	5	4.6	4.7	3	0.1	2.1	6.4	20	>6	4.6	4.6	47	0.1	2.2	4.3	20	>6
Lab043	5	4.2	4.7	4	0.1	2.1	16.2	20	3.9	4.2	4.5	85	0.5	11.1	29.4	20	1.2
Lab049	5	4.7	4.9	4	0.3	6.1	16.5	20	2.6	4.7	4.7	91	0.2	4.3	8.5	20	4.7
Total	-	-	4.7	30	0.3	6.4	-	-	-	-	4.6	650	0.3	6.5	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2023 ▼ 04 ▼ Change

Cumulative : from 2021 ▼ 01 ▼ 01 ▼ to 2023 ▼ 04 ▼ 30 ▼ Change

[TOP](#)

Reagent Kit	Reagent Code
Medicon	5