

# G6PD 定量檢驗內部品管統計結果及 SDI 分佈圖

選擇 LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

列印圖表

## 檢驗單位 : Lab009

QC Control Lot No.	Control N		Control D	
	BJ0922N		BJ0922D	
檢驗期間	Month (2023/05)	CUM (2022/12/02~2023/05/31)	Month (2023/05)	CUM (2022/12/02~2023/05/31)
累計點數 (N)	4	26	4	26
累計平均值 (Mean) U/gHb	14.2	14.3	5.8	5.9
累計標準差 (SD)	0.3	0.5	0.1	0.2
累計變異係數 (CV) %	2.1	3.5	1.7	3.4
目標值 (Target) U/gHb	16.3	16.3	5.9	5.9
總誤差 (TE) %	17.1	19.3	5.1	6.8
可容許總誤差 (TEa) %	20	20	20	20
$\sigma$	3.4	2.2	>6	5.9

誤差 (Bias)% =  $\{ [ | \text{累計平均值 (Mean)} - \text{目標值 (Target)} | ] / \text{目標值 (Target)} \} \times 100\%$   
 總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)  
 $\sigma$  (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2023 05 Change ; 累計期間 (CUM) : 自 2022 12 02 至 2023 05 31 Change


[TOP](#)

## Control N 及 Control D SDI 分佈圖

Lot No.: BJ0922N ; 檢驗期間 : 2022-12-02 至 2023-05-31 ; 目標值 (Target) : 16.3 ; 參考標準差 (SD) : 1.38 (●)

Lot No.: BJ0922D ; 檢驗期間 : 2022-12-02 至 2023-05-31 ; 目標值 (Target) : 5.9 ; 參考標準差 (SD) : 0.47 (□)

單位 : Lab009

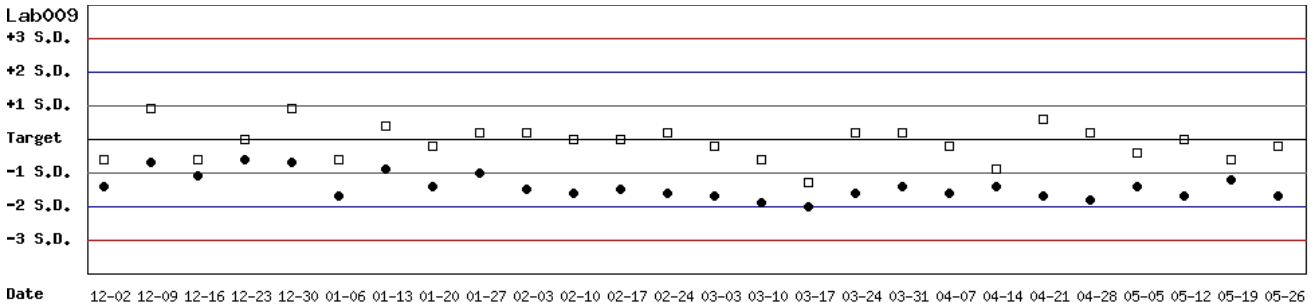


月統計 : 2023 05 Change ; 累計期間 (CUM) : 自 2022 12 02 至 2023 05 31 Change

[TOP](#)

## Control N SDI 分佈圖

Lot No.: BJ0922N ; 檢驗期間 : 2022-12-02 至 2023-05-31 ; 目標值 (Target) : 16.3 ; 參考標準差 (SD) : 1.38

單位 : Lab009

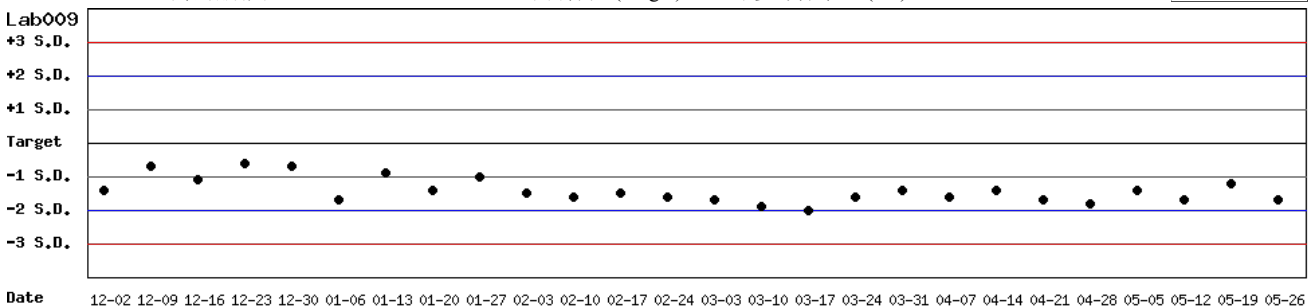


月統計 : 2023 05 Change ; 累計期間 (CUM) : 自 2022 12 02 至 2023 05 31 Change

[TOP](#)

# Control D SDI 分佈圖

Lot No.: BJ0922D ; 檢驗期間 : 2022-12-02 至 2023-05-31 ; 目標值 (Target) : 5.9 ; 參考標準差 (SD) : 0.47

單位 : Lab009

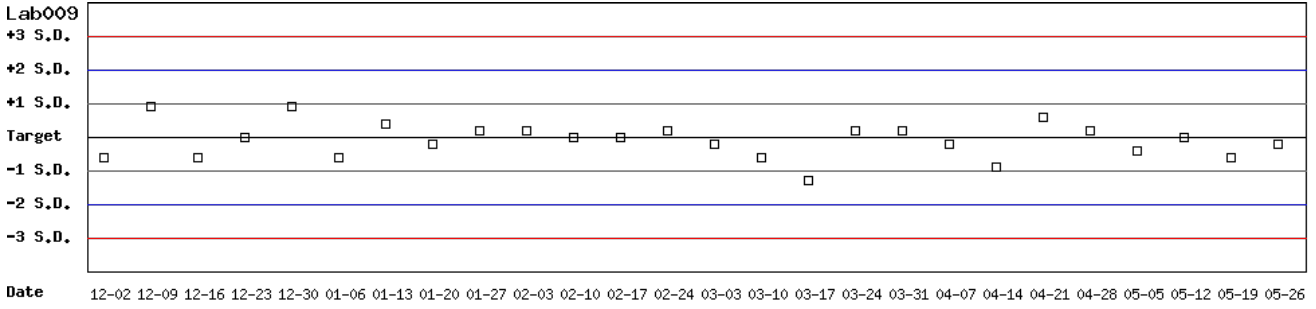

月統計 :    ; 累計期間 (CUM) : 自    至

[TOP](#)

# 各檢驗單位 G6PD 定量檢驗 內部品管統計結果 (總表一)

選擇 LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

選擇 Reagent Kit : 3 - Innovation Change

[列印總表一](#)

## 月統計表

統計月份 : 2023 05 Change

UnitID <span>↑</span>	Reagent Kit (Code) <span>↑</span>	Control N (Lot No.: BJ0922N)								Control D (Lot No.: BJ0922D)							
		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
<a href="#">Lab004</a>	3	16.2	16.1	5	0.3	1.9	4.3	20	>6	5.2	5.3	5	0.2	3.8	9.5	20	4.8
<a href="#">Lab009</a>	3	16.3	14.2	4	0.3	2.1	17.1	20	3.4	5.9	5.8	4	0.1	1.7	5.1	20	>6
<a href="#">Lab051</a>	3	15.3	14.6	7	0.4	2.7	10.1	20	5.7	5.4	7.4	7	2.9	39.2	115.4	20	-0.4
Total	-	-	15.0	16	0.8	5.3	-	-	-	-	6.3	16	2.1	33.3	-	-	-

誤差 (Bias)% = { [ | 累計平均值 (Mean) - 目標值 (Target) | ] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)% ) / 累計變異係數 (CV)

[\[TOP\]](#)

## 累計統計表

累計期間 : 自 2021 01 01 至 2023 05 31 Change

UnitID <span>↑</span>	Reagent Kit (Code) <span>↑</span>	Control N (Lot No.: BJ0922N)								Control D (Lot No.: BJ0922D)							
		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
<a href="#">Lab004</a>	3	16.2	16.2	111	0.4	2.5	4.9	20	>6	5.2	5.2	111	0.2	3.8	7.7	20	5.3
<a href="#">Lab009</a>	3	16.3	14.3	26	0.5	3.5	19.3	20	2.2	5.9	5.9	26	0.2	3.4	6.8	20	5.9
<a href="#">Lab013</a>	3	17.1	16.9	106	0.2	1.2	3.5	20	>6	5.6	5.6	106	0.2	3.6	7.1	20	5.6
<a href="#">Lab014</a>	3	16.2	15.4	113	0.9	5.8	16.6	20	2.6	5.3	5.5	113	0.2	3.6	11.0	20	4.5
<a href="#">Lab016</a>	3	13.8	14.7	72	0.4	2.7	12.0	20	5.0	4.7	4.7	72	0.2	4.3	8.5	20	4.7
<a href="#">Lab029</a>	3	15.6	14.0	34	0.3	2.1	14.5	20	4.6	5.4	5.5	34	0.2	3.6	9.1	20	5.0
<a href="#">Lab048</a>	3	15.5	15.4	21	0.4	2.6	5.8	20	>6	5.0	5.1	21	0.2	3.9	9.8	20	4.6
<a href="#">Lab051</a>	3	15.3	14.6	82	0.3	2.1	8.7	20	>6	5.4	5.8	82	0.9	15.5	38.4	20	0.8
Total	-	-	15.5	565	1.1	7.1	-	-	-	-	5.4	565	0.5	9.3	-	-	-

誤差 (Bias)% = { [ | 累計平均值 (Mean) - 目標值 (Target) | ] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)% ) / 累計變異係數 (CV)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Innovation	3

# 各檢驗單位 G6PD 定量檢驗 內部品管統計結果 (總表二)

選擇 LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

選擇 Reagent Kit : 3 - Innovation Change

列印總表二

## Control N 月統計與累計統計之比較表

Control N (Lot No.: BJ0922N)																	
Month (2023/05)										CUM (2021/01/01~2023/05/31)							
UnitID <span>↑</span>	Reagent Kit (Code) <span>↑</span>	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
<a href="#">Lab004</a>	3	16.2	16.1	5	0.3	1.9	4.3	20	>6	16.2	16.2	111	0.4	2.5	4.9	20	>6
<a href="#">Lab009</a>	3	16.3	14.2	4	0.3	2.1	17.1	20	3.4	16.3	14.3	26	0.5	3.5	19.3	20	2.2
<a href="#">Lab013</a>	3	17.1	-	0	-	-	-	20	-	17.1	16.9	106	0.2	1.2	3.5	20	>6
<a href="#">Lab014</a>	3	16.2	-	0	-	-	-	20	-	16.2	15.4	113	0.9	5.8	16.6	20	2.6
<a href="#">Lab016</a>	3	13.8	-	0	-	-	-	20	-	13.8	14.7	72	0.4	2.7	12.0	20	5.0
<a href="#">Lab029</a>	3	15.6	-	0	-	-	-	20	-	15.6	14.0	34	0.3	2.1	14.5	20	4.6
<a href="#">Lab048</a>	3	15.5	-	0	-	-	-	20	-	15.5	15.4	21	0.4	2.6	5.8	20	>6
<a href="#">Lab051</a>	3	15.3	14.6	7	0.4	2.7	10.1	20	5.7	15.3	14.6	82	0.3	2.1	8.7	20	>6
Total	-	-	15.0	16	0.8	5.3	-	-	-	-	15.5	565	1.1	7.1	-	-	-

誤差 (Bias)% = { [ | 累計平均值 (Mean) - 目標值 (Target) | ] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)% ) / 累計變異係數 (CV)

月統計 : 2023 ↓ 05 ↓ Change

累計期間 (CUM) : 自 2021 ↓ 01 ↓ 01 ↓ 至 2023 ↓ 05 ↓ 31 ↓ Change

[TOP](#)

## Control D 月統計與累計統計之比較表

Control D (Lot No.: BJ0922D)																	
Month (2023/05)										CUM (2021/01/01~2023/05/31)							
UnitID <span>↑</span>	Reagent Kit (Code) <span>↑</span>	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
<a href="#">Lab004</a>	3	5.2	5.3	5	0.2	3.8	9.5	20	4.8	5.2	5.2	111	0.2	3.8	7.7	20	5.3
<a href="#">Lab009</a>	3	5.9	5.8	4	0.1	1.7	5.1	20	>6	5.9	5.9	26	0.2	3.4	6.8	20	5.9
<a href="#">Lab013</a>	3	5.6	-	0	-	-	-	20	-	5.6	5.6	106	0.2	3.6	7.1	20	5.6
<a href="#">Lab014</a>	3	5.3	-	0	-	-	-	20	-	5.3	5.5	113	0.2	3.6	11.0	20	4.5
<a href="#">Lab016</a>	3	4.7	-	0	-	-	-	20	-	4.7	4.7	72	0.2	4.3	8.5	20	4.7
<a href="#">Lab029</a>	3	5.4	-	0	-	-	-	20	-	5.4	5.5	34	0.2	3.6	9.1	20	5.0
<a href="#">Lab048</a>	3	5.0	-	0	-	-	-	20	-	5.0	5.1	21	0.2	3.9	9.8	20	4.6
<a href="#">Lab051</a>	3	5.4	7.4	7	2.9	39.2	115.4	20	-0.4	5.4	5.8	82	0.9	15.5	38.4	20	0.8
Total	-	-	6.3	16	2.1	33.3	-	-	-	-	5.4	565	0.5	9.3	-	-	-

誤差 (Bias)% = { [ | 累計平均值 (Mean) - 目標值 (Target) | ] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)% ) / 累計變異係數 (CV)

月統計 : 2023 ↓ 05 ↓ Change

累計期間 (CUM) : 自 2021 ↓ 01 ↓ 01 ↓ 至 2023 ↓ 05 ↓ 31 ↓ Change

[TOP](#)

Reagent Kit	Reagent Code
Innovation	3