

QC Chart of Internal Quality Control (IQC)for G6PD Quantitative Test

Select LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

Print Table

Lab037

	Control N		Control D	
QC Control Lot No.	BA1116N		BA1116D	
Duration of the Analyzing	Month (2025/01)	CUM (2023/11/24~2025/01/31)	Month (2025/01)	CUM (2023/11/24~2025/01/31)
Runs (N)	1	26	1	26
Mean (U/gHb)	13.2	13.7	4.5	4.4
SD	-	0.4	-	0.2
CV (%)	-	2.9	-	4.5
Target Value (U/gHb)	12.8	12.8	4.2	4.2
Total Error (%)	-	12.9	-	13.9
TEa (%)	20	20	20	20
σ	-	4.5	-	3.4

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2025 01 Change ; Cumulative : from 2023 11 24 to 2025 01 31 Change


[\[TOP\]](#)

Control N and Control D SDI QC Chart

Lot No.: BA1116N ; Duration : 2023-11-24 to 2025-01-31 ; Target : 12.8 ; SD : 1.28 (●)

Lot No.: BA1116D ; Duration : 2023-11-24 to 2025-01-31 ; Target : 4.2 ; SD : 0.42 (□)

Lab037

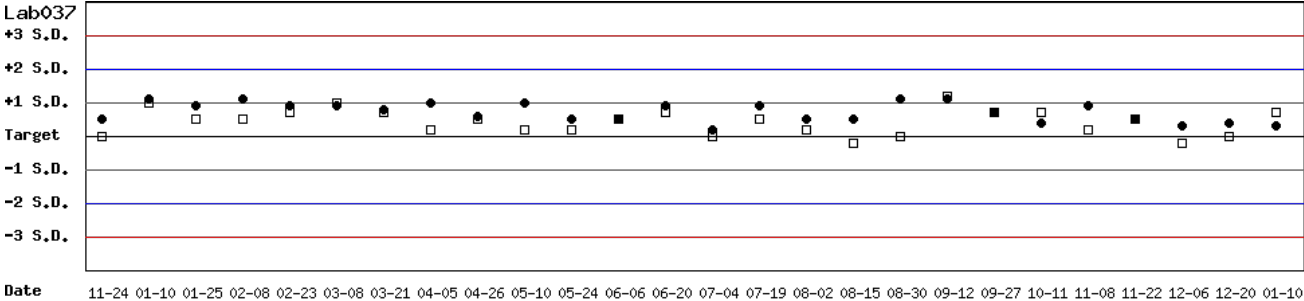


Month : 2025 01 Change ; Cumulative : from 2023 11 24 to 2025 01 31 Change

[\[TOP\]](#)

Control N SDI QC Chart

Lot No.: BA1116N ; Duration : 2023-11-24 to 2025-01-31 ; Target : 12.8 ; SD : 1.28

Lab037

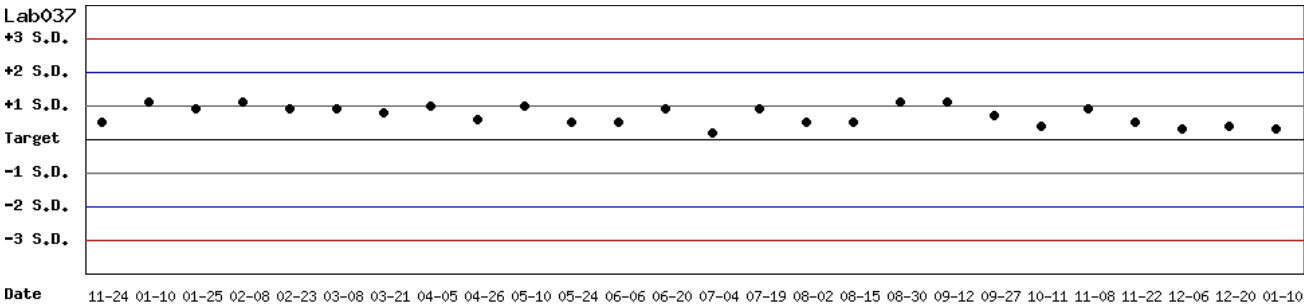

Month : 2025 01 Change ; Cumulative : from 2023 11 24 to 2025 01 31 Change

[\[TOP\]](#)

Control D SDI QC Chart

Lab037

Lot No.: BA1116D ; Duration : 2023-11-24 to 2025-01-31 ; Target : 4.2 ; SD : 0.42

Lab037

+3 S.D.

+2 S.D.

+1 S.D.

Target

-1 S.D.

-2 S.D.

-3 S.D.

Date

11-24 01-10 01-25 02-08 02-23 03-08 03-21 04-05 04-26 05-10 05-24 06-06 06-20 07-04 07-19 08-02 08-15 08-30 09-12 09-27 10-11 11-08 11-22 12-06 12-20 01-10

Month : 2025 01 Change ; Cumulative : from 2023 11 24 to 2025 01 31 Change

[TOP]

Peer Group Statistics (Table 1)

Select LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

[Print Table 1](#)

Monthly

Month : 2025 01 Change

		Control N (Lot No.: BA1116N)									Control D (Lot No.: BA1116D)							
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	13.6	14.1	9	0.1	0.7	5.1	20	>6		4.3	4.4	9	0.1	2.3	6.9	20	>6
Lab033	5	15.8	15.3	5	0.2	1.3	5.8	20	>6		4.7	4.6	5	0.2	4.3	10.8	20	4.2
Lab037	5	12.8	13.2	1	-	-	-	20	-		4.2	4.5	1	-	-	-	20	-
Lab052	5	15.3	15.3	1	-	-	-	20	-		4.6	4.7	1	-	-	-	20	-
Total	-	-	14.5	16	0.7	4.8	-	-	-		-	4.5	16	0.2	4.4	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Cumulative

Cumulative : from 2022 06 01 to 2025 01 31 Change

		Control N (Lot No.: BA1116N)									Control D (Lot No.: BA1116D)							
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	13.6	13.8	139	0.4	2.9	7.3	20	>6		4.3	4.4	139	0.1	2.3	6.9	20	>6
Lab027	5	12.8	14.9	74	0.9	6.0	28.5	20	0.6		4.2	4.9	74	0.4	8.2	33.0	20	0.4
Lab032	5	14.9	13.7	22	0.4	2.9	13.9	20	4.1		5.1	4.2	22	0.2	4.8	27.2	20	0.5
Lab033	5	15.8	15.7	79	0.5	3.2	7.0	20	>6		4.7	4.8	79	0.2	4.2	10.5	20	4.3
Lab037	5	12.8	13.8	33	0.4	2.9	13.6	20	4.2		4.2	4.4	33	0.2	4.5	13.9	20	3.4
Lab052	5	15.3	15.4	9	0.2	1.3	3.3	20	>6		4.6	4.7	9	0.1	2.1	6.4	20	>6
Total	-	-	14.5	360	1.0	6.9	-	-	-		-	4.6	360	0.3	6.5	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Medicon	5

Peer Group Statistics (Table 2)

Select LotNo :
BA1116N (2022-04-26 ~ 2100-12-31)
Change

Select Reagent Kit :
5 - Medicon
Change

Print Table 2

Control N Month vs. Cumulative

		Control N (Lot No.: BA1116N)																	
		Month (2025/01)										CUM (2022/06/01~2025/01/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	
Lab018	5	13.6	14.1	9	0.1	0.7	5.1	20	>6		13.6	13.8	139	0.4	2.9	7.3	20	>6	
Lab027	5	12.8	-	0	-	-	-	20	-		12.8	14.9	74	0.9	6.0	28.5	20	0.6	
Lab032	5	14.9	-	0	-	-	-	20	-		14.9	13.7	22	0.4	2.9	13.9	20	4.1	
Lab033	5	15.8	15.3	5	0.2	1.3	5.8	20	>6		15.8	15.7	79	0.5	3.2	7.0	20	>6	
Lab037	5	12.8	13.2	1	-	-	-	20	-		12.8	13.8	33	0.4	2.9	13.6	20	4.2	
Lab052	5	15.3	15.3	1	-	-	-	20	-		15.3	15.4	9	0.2	1.3	3.3	20	>6	
Total	-	-	14.5	16	0.7	4.8	-	-	-		-	14.5	360	1.0	6.9	-	-	-	

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month :
2025
01
Change
Cumulative : from
2022
06
01
to
2025
01
31
Change

[TOP](#)

Control D Month vs. Cumulative

		Control D (Lot No.: BA1116D)																
		Month (2025/01)									CUM (2022/06/01~2025/01/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	4.3	4.4	9	0.1	2.3	6.9	20	>6		4.3	4.4	139	0.1	2.3	6.9	20	>6
Lab027	5	4.2	-	0	-	-	-	20	-		4.2	4.9	74	0.4	8.2	33.0	20	0.4
Lab032	5	5.1	-	0	-	-	-	20	-		5.1	4.2	22	0.2	4.8	27.2	20	0.5
Lab033	5	4.7	4.6	5	0.2	4.3	10.8	20	4.2		4.7	4.8	79	0.2	4.2	10.5	20	4.3
Lab037	5	4.2	4.5	1	-	-	-	20	-		4.2	4.4	33	0.2	4.5	13.9	20	3.4
Lab052	5	4.6	4.7	1	-	-	-	20	-		4.6	4.7	9	0.1	2.1	6.4	20	>6
Total	-	-	4.5	16	0.2	4.4	-	-	-		-	4.6	360	0.3	6.5	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month :
2025
01
Change
Cumulative : from
2022
06
01
to
2025
01
31
Change

[TOP](#)

Reagent Kit	Reagent Code
Medicon	5