

G6PD 定量檢驗內部品管統計結果及 SDI 分佈圖

選擇 LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

列印圖表

檢驗單位 : Lab006

	Control N		Control D	
QC Control Lot No.	BA1116N		BA1116D	
檢驗期間	Month (2025/02)	CUM (2024/11/27~2025/02/28)	Month (2025/02)	CUM (2024/11/27~2025/02/28)
累計點數 (N)	8	26	8	26
累計平均值 (Mean) U/gHb	14.6	14.0	5.0	4.7
累計標準差 (SD)	0.4	0.6	0.2	0.3
累計變異係數 (CV) %	2.7	4.3	4.0	6.4
目標值 (Target) U/gHb	14.5	14.5	4.9	4.9
總誤差 (TE) %	6.2	12.0	10.0	16.8
可容許總誤差 (TEa) %	20	20	20	20
σ	>6	3.8	4.5	2.5

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2025 02 Change ; 累計期間 (CUM) : 自 2024 11 27 至 2025 02 28 Change


TOP

Control N 及 Control D SDI
分佈圖

Lot No.: BA1116N; 檢驗期間 : 2024-11-27 至 2025-02-28; 目標值 (Target): 14.5; 參考標準差 (SD): 1.00 (●)

Lot No.: BA1116D; 檢驗期間 : 2024-11-27 至 2025-02-28; 目標值 (Target): 4.9; 參考標準差 (SD): 0.20 (□)

單位 : Lab006

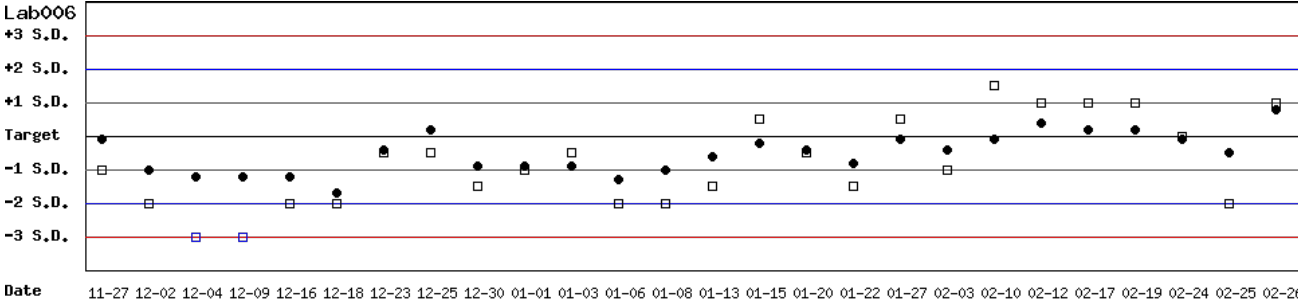


月統計 : 2025 02 Change ; 累計期間 (CUM) : 自 2024 11 27 至 2025 02 28 Change

TOP

Control N SDI 分佈圖

Lot No.: BA1116N; 檢驗期間 : 2024-11-27 至 2025-02-28; 目標值 (Target): 14.5; 參考標準差 (SD): 1.00

單位 : Lab006

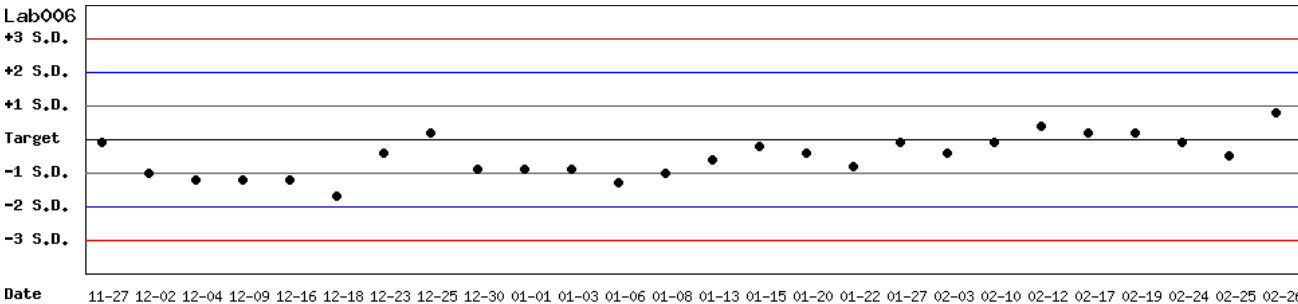


月統計 : 2025 02 Change ; 累計期間 (CUM) : 自 2024 11 27 至 2025 02 28 Change

TOP

Control D SDI 分佈圖

Lot No.: BA1116D ; 檢驗期間 : 2024-11-27 至 2025-02-28 ; 目標值 (Target) : 4.9 ; 參考標準差 (SD) : 0.20

單位 : Lab006

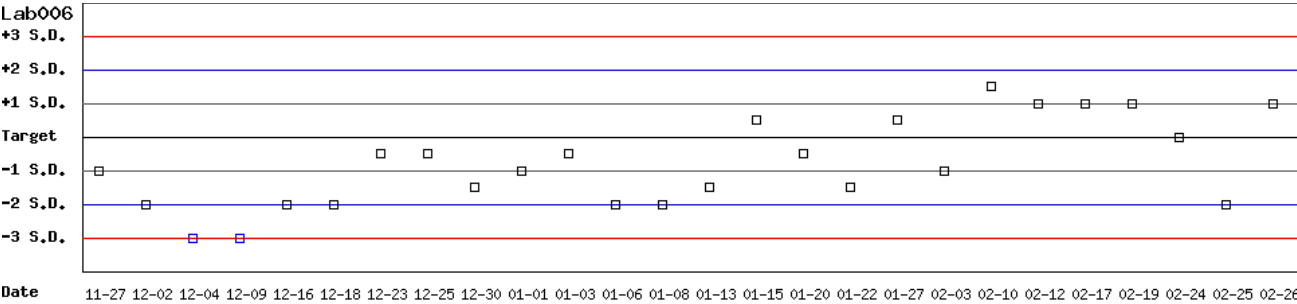

Date 11-27 12-02 12-04 12-09 12-16 12-18 12-23 12-25 12-30 01-01 01-03 01-06 01-08 01-13 01-15 01-20 01-22 01-27 02-03 02-10 02-12 02-17 02-19 02-24 02-25 02-26

月統計 : 2025 02 Change ; 累計期間 (CUM) : 自 2024 11 27 至 2025 02 28 Change

[\[TOP\]](#)

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表一）

選擇 LotNo： BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit： 4 - Lanner Change

列印總表一

月統計表

統計月份： 2025 02 Change

		Control N (Lot No.: BA1116N)									Control D (Lot No.: BA1116D)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	14.5	14.6	8	0.4	2.7	6.2	20	>6		4.9	5.0	8	0.2	4.0	10.0	20	4.5
Lab007	4	16.8	17.1	8	0.4	2.3	6.5	20	>6		5.1	5.1	8	0.1	2.0	3.9	20	>6
Lab008	4	14.2	13.9	8	0.4	2.9	7.9	20	>6		4.5	4.7	8	0.2	4.3	13.0	20	3.6
Lab010	4	14.9	14.9	1	-	-	-	20	-		4.8	5.2	1	-	-	-	20	-
Total	-	-	15.2	25	1.4	9.2	-	-	-		-	4.9	25	0.3	6.1	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[\[TOP\]](#)

累計統計表

累計期間：自 2022 06 01 至 2025 02 28 Change

		Control N (Lot No.: BA1116N)										Control D (Lot No.: BA1116D)							
UnitID 	Reagent Kit (Code) 	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	
Lab006	4	14.5	14.5	130	0.5	3.4	6.9	20	5.9		4.9	4.8	130	0.3	6.3	14.5	20	2.9	
Lab007	4	16.8	15.5	215	1.0	6.5	20.6	20	1.9		5.1	5.0	215	0.4	8.0	18.0	20	2.3	
Lab008	4	14.2	14.1	121	0.6	4.3	9.2	20	4.5		4.5	4.7	121	0.3	6.4	17.2	20	2.4	
Lab010	4	14.9	15.0	56	0.4	2.7	6.0	20	>6		4.8	4.8	56	0.5	10.4	20.8	20	1.9	
Lab011	4	15.2	14.6	21	0.3	2.1	8.1	20	>6		5.4	3.9	21	0.2	5.1	38.0	20	-1.5	
Lab012	4	16.0	15.0	166	0.8	5.3	16.9	20	2.6		5.1	5.1	166	0.4	7.8	15.7	20	2.6	
Total	-	-	14.9	709	0.9	6.0	-	-	-		-	4.9	709	0.4	8.2	-	-	-	

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Lanner	4

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表二）

選擇 LotNo：BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit：4 - Lanner Change

列印總表二

Control N 月統計與累計統計之比較表

		Control N (Lot No.: BA1116N)																
		Month (2025/02)									CUM (2022/06/01~2025/02/28)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	14.5	14.6	8	0.4	2.7	6.2	20	>6		14.5	14.5	130	0.5	3.4	6.9	20	5.9
Lab007	4	16.8	17.1	8	0.4	2.3	6.5	20	>6		16.8	15.5	215	1.0	6.5	20.6	20	1.9
Lab008	4	14.2	13.9	8	0.4	2.9	7.9	20	>6		14.2	14.1	121	0.6	4.3	9.2	20	4.5
Lab010	4	14.9	14.9	1	-	-	-	20	-		14.9	15.0	56	0.4	2.7	6.0	20	>6
Lab011	4	15.2	-	0	-	-	-	20	-		15.2	14.6	21	0.3	2.1	8.1	20	>6
Lab012	4	16.0	-	0	-	-	-	20	-		16.0	15.0	166	0.8	5.3	16.9	20	2.6
Total	-	-	15.2	25	1.4	9.2	-	-	-		-	14.9	709	0.9	6.0	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計：2025 02 Change 累計期間 (CUM)：自 2022 06 01 至 2025 02 28 Change

[TOP](#)

Control D 月統計與累計統計之比較表

		Control D (Lot No.: BA1116D)																
		Month (2025/02)									CUM (2022/06/01~2025/02/28)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	4.9	5.0	8	0.2	4.0	10.0	20	4.5		4.9	4.8	130	0.3	6.3	14.5	20	2.9
Lab007	4	5.1	5.1	8	0.1	2.0	3.9	20	>6		5.1	5.0	215	0.4	8.0	18.0	20	2.3
Lab008	4	4.5	4.7	8	0.2	4.3	13.0	20	3.6		4.5	4.7	121	0.3	6.4	17.2	20	2.4
Lab010	4	4.8	5.2	1	-	-	-	20	-		4.8	4.8	56	0.5	10.4	20.8	20	1.9
Lab011	4	5.4	-	0	-	-	-	20	-		5.4	3.9	21	0.2	5.1	38.0	20	-1.5
Lab012	4	5.1	-	0	-	-	-	20	-		5.1	5.1	166	0.4	7.8	15.7	20	2.6
Total	-	-	4.9	25	0.3	6.1	-	-	-		-	4.9	709	0.4	8.2	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計：2025 02 Change 累計期間 (CUM)：自 2022 06 01 至 2025 02 28 Change

[TOP](#)

Reagent Kit	Reagent Code
Lanner	4