

G6PD 定量檢驗內部品管統計結果及 SDI 分佈圖

選擇 LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

列印圖表

檢驗單位 : Lab008

	Control N		Control D	
QC Control Lot No.	BA1116N		BA1116D	
檢驗期間	Month (2025/02)	CUM (2024/11/27~2025/02/28)	Month (2025/02)	CUM (2024/11/27~2025/02/28)
累計點數 (N)	8	26	8	26
累計平均值 (Mean) U/gHb	13.9	13.4	4.7	4.5
累計標準差 (SD)	0.4	0.6	0.2	0.3
累計變異係數 (CV) %	2.9	4.5	4.3	6.7
目標值 (Target) U/gHb	14.2	14.2	4.5	4.5
總誤差 (TE) %	7.9	14.6	13.0	13.3
可容許總誤差 (TEa) %	20	20	20	20
σ	>6	3.2	3.6	3.0

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2025 02 Change ; 累計期間 (CUM) : 自 2024 11 27 至 2025 02 28 Change


TOP

Control N 及 Control D SDI
分佈圖

Lot No.: BA1116N; 檢驗期間 : 2024-11-27 至 2025-02-28; 目標值 (Target): 14.2; 參考標準差 (SD): 1.28 (●)

Lot No.: BA1116D; 檢驗期間 : 2024-11-27 至 2025-02-28; 目標值 (Target): 4.5; 參考標準差 (SD): 0.42 (□)

單位 : Lab008

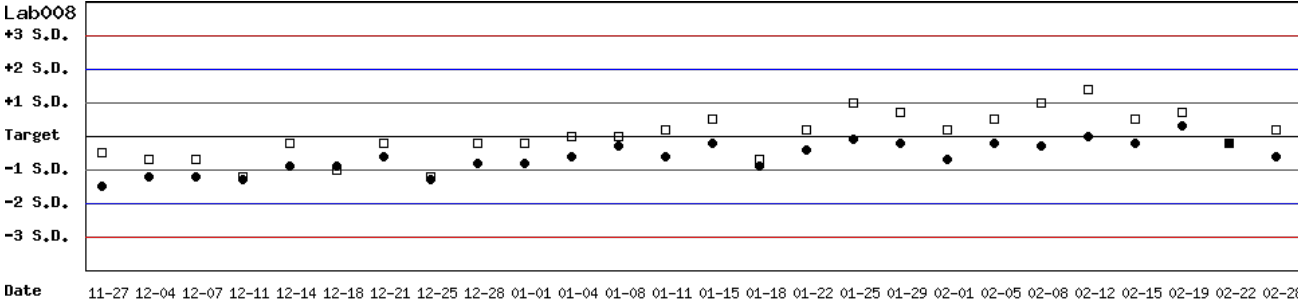


月統計 : 2025 02 Change ; 累計期間 (CUM) : 自 2024 11 27 至 2025 02 28 Change

TOP

Control N SDI 分佈圖

Lot No.: BA1116N; 檢驗期間 : 2024-11-27 至 2025-02-28; 目標值 (Target): 14.2; 參考標準差 (SD): 1.28

單位 : Lab008

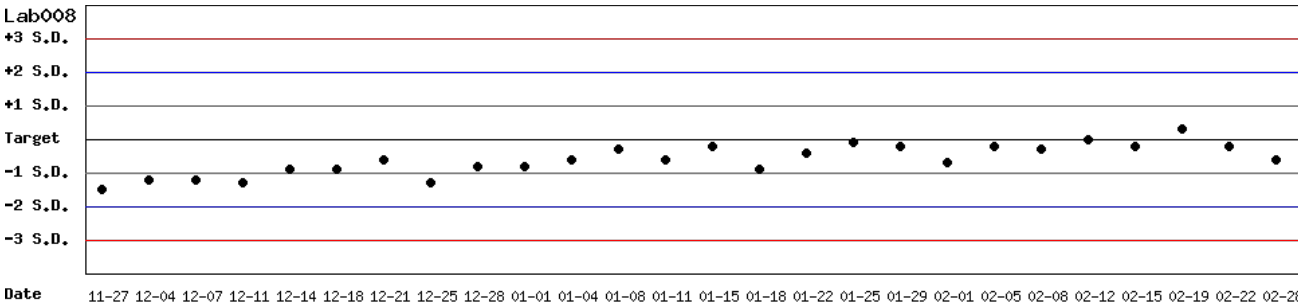


月統計 : 2025 02 Change ; 累計期間 (CUM) : 自 2024 11 27 至 2025 02 28 Change

TOP

Control D SDI 分佈圖

單位: Lab008

Lot No.: BA1116D; 檢驗期間: 2024-11-27 至 2025-02-28; 目標值 (Target): 4.5; 參考標準差 (SD): 0.42

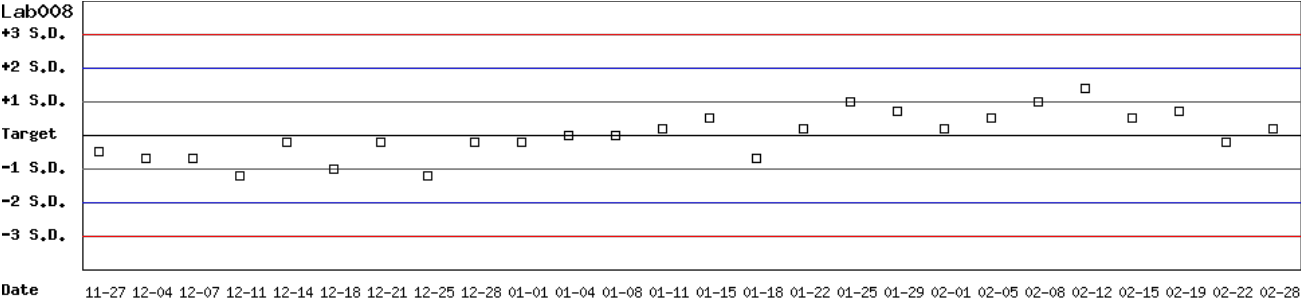

Date 11-27 12-04 12-07 12-11 12-14 12-18 12-21 12-25 12-28 01-01 01-04 01-08 01-11 01-15 01-18 01-22 01-25 01-29 02-01 02-05 02-08 02-12 02-15 02-19 02-22 02-28

月統計: 2025 02 Change ; 累計期間 (CUM): 自 2024 11 27 至 2025 02 28 Change

[\[TOP\]](#)

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表一）

選擇 LotNo： BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit： 4 - Lanner Change

列印總表一

月統計表

統計月份： 2025 02 Change

		Control N (Lot No.: BA1116N)									Control D (Lot No.: BA1116D)								
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	
Lab006	4	14.5	14.6	8	0.4	2.7	6.2	20	>6		4.9	5.0	8	0.2	4.0	10.0	20	4.5	
Lab007	4	16.8	17.1	8	0.4	2.3	6.5	20	>6		5.1	5.1	8	0.1	2.0	3.9	20	>6	
Lab008	4	14.2	13.9	8	0.4	2.9	7.9	20	>6		4.5	4.7	8	0.2	4.3	13.0	20	3.6	
Lab010	4	14.9	14.9	1	-	-	-	20	-		4.8	5.2	1	-	-	-	20	-	
Total	-	-	15.2	25	1.4	9.2	-	-	-		-	4.9	25	0.3	6.1	-	-	-	

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[TOP]

累計統計表

累計期間：自 2022 06 01 至 2025 02 28 Change

		Control N (Lot No.: BA1116N)									Control D (Lot No.: BA1116D)								
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	
Lab006	4	14.5	14.5	130	0.5	3.4	6.9	20	5.9		4.9	4.8	130	0.3	6.3	14.5	20	2.9	
Lab007	4	16.8	15.5	215	1.0	6.5	20.6	20	1.9		5.1	5.0	215	0.4	8.0	18.0	20	2.3	
Lab008	4	14.2	14.1	121	0.6	4.3	9.2	20	4.5		4.5	4.7	121	0.3	6.4	17.2	20	2.4	
Lab010	4	14.9	15.0	56	0.4	2.7	6.0	20	>6		4.8	4.8	56	0.5	10.4	20.8	20	1.9	
Lab011	4	15.2	14.6	21	0.3	2.1	8.1	20	>6		5.4	3.9	21	0.2	5.1	38.0	20	-1.5	
Lab012	4	16.0	15.0	166	0.8	5.3	16.9	20	2.6		5.1	5.1	166	0.4	7.8	15.7	20	2.6	
Total	-	-	14.9	709	0.9	6.0	-	-	-		-	4.9	709	0.4	8.2	-	-	-	

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[TOP]

Reagent Kit	Reagent Code
Lanner	4

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表二）

選擇 LotNo： BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit： 4 - Lanner Change

列印總表二

Control N 月統計與累計統計之比較表

Control N (Lot No.: BA1116N)																	
Month (2025/02)										CUM (2022/06/01~2025/02/28)							
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	14.5	14.6	8	0.4	2.7	6.2	20	>6	14.5	14.5	130	0.5	3.4	6.9	20	5.9
Lab007	4	16.8	17.1	8	0.4	2.3	6.5	20	>6	16.8	15.5	215	1.0	6.5	20.6	20	1.9
Lab008	4	14.2	13.9	8	0.4	2.9	7.9	20	>6	14.2	14.1	121	0.6	4.3	9.2	20	4.5
Lab010	4	14.9	14.9	1	-	-	-	20	-	14.9	15.0	56	0.4	2.7	6.0	20	>6
Lab011	4	15.2	-	0	-	-	-	20	-	15.2	14.6	21	0.3	2.1	8.1	20	>6
Lab012	4	16.0	-	0	-	-	-	20	-	16.0	15.0	166	0.8	5.3	16.9	20	2.6
Total	-	-	15.2	25	1.4	9.2	-	-	-	-	14.9	709	0.9	6.0	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計： 2025 02 Change 累計期間 (CUM)：自 2022 06 01 至 2025 02 28 Change

[\[TOP\]](#)

Control D 月統計與累計統計之比較表

Control D (Lot No.: BA1116D)																	
Month (2025/02)										CUM (2022/06/01~2025/02/28)							
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	4.9	5.0	8	0.2	4.0	10.0	20	4.5	4.9	4.8	130	0.3	6.3	14.5	20	2.9
Lab007	4	5.1	5.1	8	0.1	2.0	3.9	20	>6	5.1	5.0	215	0.4	8.0	18.0	20	2.3
Lab008	4	4.5	4.7	8	0.2	4.3	13.0	20	3.6	4.5	4.7	121	0.3	6.4	17.2	20	2.4
Lab010	4	4.8	5.2	1	-	-	-	20	-	4.8	4.8	56	0.5	10.4	20.8	20	1.9
Lab011	4	5.4	-	0	-	-	-	20	-	5.4	3.9	21	0.2	5.1	38.0	20	-1.5
Lab012	4	5.1	-	0	-	-	-	20	-	5.1	5.1	166	0.4	7.8	15.7	20	2.6
Total	-	-	4.9	25	0.3	6.1	-	-	-	-	4.9	709	0.4	8.2	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計： 2025 02 Change 累計期間 (CUM)：自 2022 06 01 至 2025 02 28 Change

[\[TOP\]](#)

Reagent Kit	Reagent Code
Lanner	4