

G6PD 定量檢驗內部品管統計結果及 SDI 分佈圖

選擇 LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

列印圖表

檢驗單位 : Lab008

	Control N		Control D	
QC Control Lot No.	BA1116N		BA1116D	
檢驗期間	Month (2025/08)	CUM (2025/03/26~2025/08/29)	Month (2025/08)	CUM (2025/03/26~2025/08/29)
累計點數 (N)	7	26	7	26
累計平均值 (Mean) U/gHb	13.8	14.2	4.4	4.6
累計標準差 (SD)	1.2	0.7	0.5	0.3
累計變異係數 (CV) %	8.7	4.9	11.4	6.5
目標值 (Target) U/gHb	14.2	14.2	4.5	4.5
總誤差 (TE) %	20.2	9.9	24.9	15.3
可容許總誤差 (TEa) %	20	20	20	20
σ	2.0	4.1	1.6	2.7

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2025 08 Change ; 累計期間 (CUM) : 自 2025 03 26 至 2025 08 29 Change


TOP

Control N 及 Control D SDI 分佈圖

Lot No.: BA1116N ; 檢驗期間 : 2025-03-26 至 2025-08-29 ; 目標值 (Target) : 14.2 ; 參考標準差 (SD) : 1.28 (●)

Lot No.: BA1116D ; 檢驗期間 : 2025-03-26 至 2025-08-29 ; 目標值 (Target) : 4.5 ; 參考標準差 (SD) : 0.42 (□)

單位 : Lab008

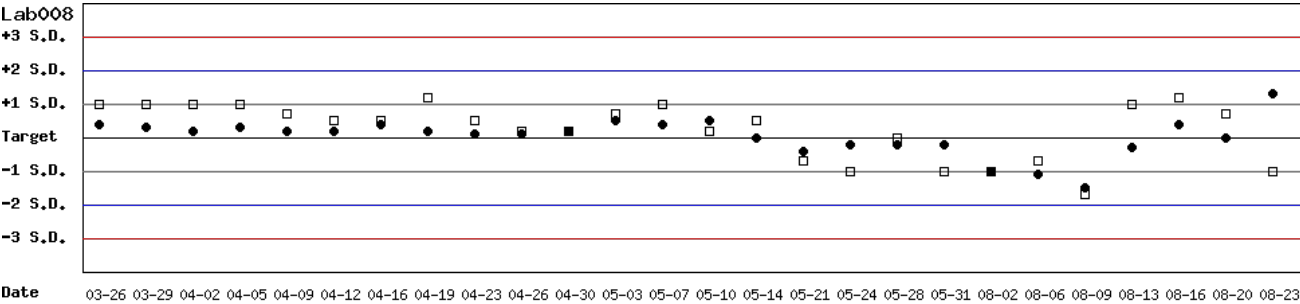


月統計 : 2025 08 Change ; 累計期間 (CUM) : 自 2025 03 26 至 2025 08 29 Change

TOP

Control N SDI 分佈圖

Lot No.: BA1116N ; 檢驗期間 : 2025-03-26 至 2025-08-29 ; 目標值 (Target) : 14.2 ; 參考標準差 (SD) : 1.28

單位 : Lab008

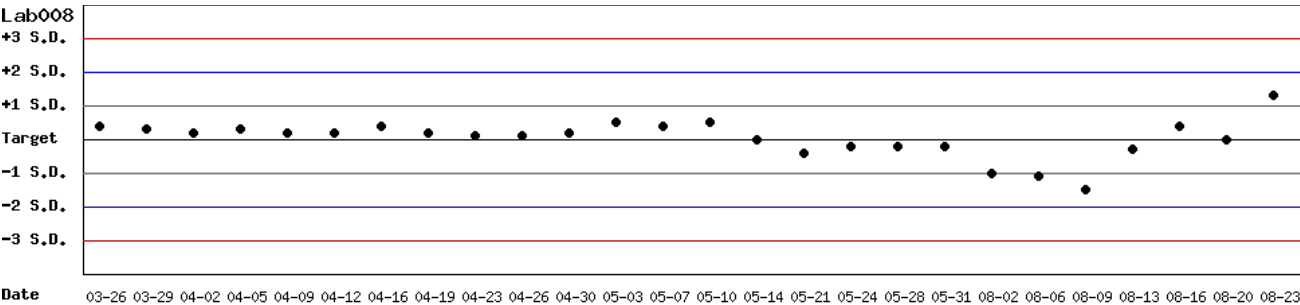


月統計 : 2025 08 Change ; 累計期間 (CUM) : 自 2025 03 26 至 2025 08 29 Change

TOP

Control D SDI 分佈圖

Lot No.: BA1116D; 檢驗期間: 2025-03-26 至 2025-08-29; 目標值 (Target): 4.5; 參考標準差 (SD): 0.42

單位: Lab008

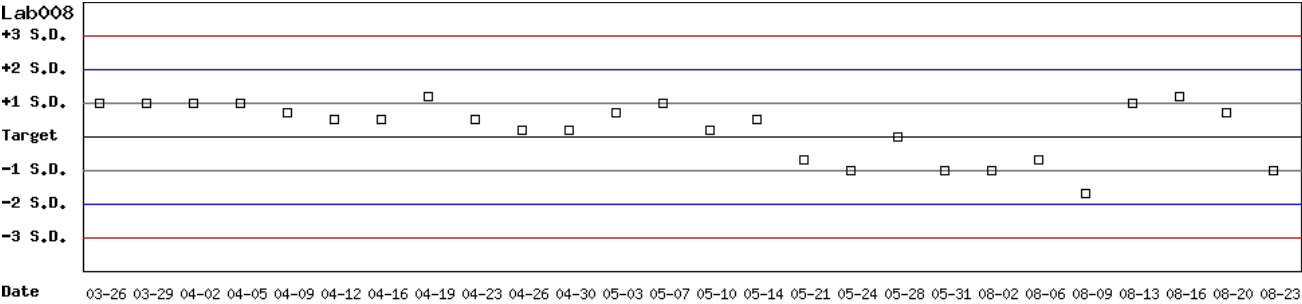

月統計: 2025 08 Change ; 累計期間 (CUM): 自 2025 03 26 至 2025 08 29 Change

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表一）

選擇 LotNo： BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit： 4 - Lanner Change

列印總表一

月統計表

統計月份： 2025 08 Change

		Control N (Lot No.: BA1116N)								Control D (Lot No.: BA1116D)							
UnitID	Reagent Kit (Code)	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab007	4	16.3	17.6	8	0.4	2.3	12.5	20	5.2	5.1	5.3	8	0.1	1.9	7.7	20	>6
Lab008	4	14.2	13.8	7	1.2	8.7	20.2	20	2.0	4.5	4.4	7	0.5	11.4	24.9	20	1.6
Lab010	4	14.9	15.7	2	-	-	-	20	-	4.8	4.8	2	-	-	-	20	-
Total	-	-	15.8	17	2.0	12.7	-	-	-	-	4.9	17	0.5	10.2	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[\[TOP\]](#)

累計統計表

累計期間：自 2022 06 01 至 2025 08 31 Change

		Control N (Lot No.: BA1116N)								Control D (Lot No.: BA1116D)							
UnitID	Reagent Kit (Code)	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	14.5	14.5	130	0.5	3.4	6.9	20	5.9	4.9	4.8	130	0.3	6.3	14.5	20	2.9
Lab007	4	16.3	15.7	265	1.0	6.4	16.4	20	2.5	5.1	5.0	265	0.4	8.0	18.0	20	2.3
Lab008	4	14.2	14.1	154	0.6	4.3	9.2	20	4.5	4.5	4.7	154	0.3	6.4	17.2	20	2.4
Lab010	4	14.9	15.0	70	0.4	2.7	6.0	20	>6	4.8	4.8	70	0.5	10.4	20.8	20	1.9
Lab011	4	15.2	14.6	21	0.3	2.1	8.1	20	>6	5.4	3.9	21	0.2	5.1	38.0	20	-1.5
Lab012	4	16.0	15.0	166	0.8	5.3	16.9	20	2.6	5.1	5.1	166	0.4	7.8	15.7	20	2.6
Total	-	-	15.0	806	1.0	6.7	-	-	-	-	4.9	806	0.4	8.2	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Lanner	4

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表二）

選擇 LotNo： BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit： 4 - Lanner Change

列印總表二

Control N 月統計與累計統計之比較表

Control N (Lot No.: BA1116N)																	
Month (2025/08)										CUM (2022/06/01~2025/08/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	14.5	-	0	-	-	-	20	-	14.5	14.5	130	0.5	3.4	6.9	20	5.9
Lab007	4	16.3	17.6	8	0.4	2.3	12.5	20	5.2	16.3	15.7	265	1.0	6.4	16.4	20	2.5
Lab008	4	14.2	13.8	7	1.2	8.7	20.2	20	2.0	14.2	14.1	154	0.6	4.3	9.2	20	4.5
Lab010	4	14.9	15.7	2	-	-	-	20	-	14.9	15.0	70	0.4	2.7	6.0	20	>6
Lab011	4	15.2	-	0	-	-	-	20	-	15.2	14.6	21	0.3	2.1	8.1	20	>6
Lab012	4	16.0	-	0	-	-	-	20	-	16.0	15.0	166	0.8	5.3	16.9	20	2.6
Total	-	-	15.8	17	2.0	12.7	-	-	-	-	15.0	806	1.0	6.7	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計： 2025 08 Change 累計期間 (CUM)：自 2022 06 01 至 2025 08 31 Change

[TOP](#)

Control D 月統計與累計統計之比較表

Control D (Lot No.: BA1116D)																	
Month (2025/08)										CUM (2022/06/01~2025/08/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab006	4	4.9	-	0	-	-	-	20	-	4.9	4.8	130	0.3	6.3	14.5	20	2.9
Lab007	4	5.1	5.3	8	0.1	1.9	7.7	20	>6	5.1	5.0	265	0.4	8.0	18.0	20	2.3
Lab008	4	4.5	4.4	7	0.5	11.4	24.9	20	1.6	4.5	4.7	154	0.3	6.4	17.2	20	2.4
Lab010	4	4.8	4.8	2	-	-	-	20	-	4.8	4.8	70	0.5	10.4	20.8	20	1.9
Lab011	4	5.4	-	0	-	-	-	20	-	5.4	3.9	21	0.2	5.1	38.0	20	-1.5
Lab012	4	5.1	-	0	-	-	-	20	-	5.1	5.1	166	0.4	7.8	15.7	20	2.6
Total	-	-	4.9	17	0.5	10.2	-	-	-	-	4.9	806	0.4	8.2	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計： 2025 08 Change 累計期間 (CUM)：自 2022 06 01 至 2025 08 31 Change

[TOP](#)

Reagent Kit	Reagent Code
Lanner	4