

G6PD 定量檢驗內部品管統計結果及 SDI 分佈圖

選擇 LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

列印圖表

檢驗單位 : Lab009

	Control N		Control D	
QC Control Lot No.	BA1116N		BA1116D	
檢驗期間	Month (2025/08)	CUM (2025/03/07~2025/08/29)	Month (2025/08)	CUM (2025/03/07~2025/08/29)
累計點數 (N)	5	26	5	26
累計平均值 (Mean) U/gHb	14.4	14.3	5.2	5.0
累計標準差 (SD)	0.1	0.4	0.1	0.2
累計變異係數 (CV) %	0.7	2.8	1.9	4.0
目標值 (Target) U/gHb	12.8	12.8	4.2	4.2
總誤差 (TE) %	13.9	17.3	27.7	27.0
可容許總誤差 (TEa) %	20	20	20	20
σ	>6	3.0	-2.0	0.2

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2025 08 Change ; 累計期間 (CUM) : 自 2025 03 07 至 2025 08 29 Change


TOP

Control N 及 Control D SDI 分佈圖

Lot No.: BA1116N ; 檢驗期間 : 2025-03-07 至 2025-08-29 ; 目標值 (Target) : 12.8 ; 參考標準差 (SD) : 1.28 (●)

Lot No.: BA1116D ; 檢驗期間 : 2025-03-07 至 2025-08-29 ; 目標值 (Target) : 4.2 ; 參考標準差 (SD) : 0.42 (□)

單位 : Lab009

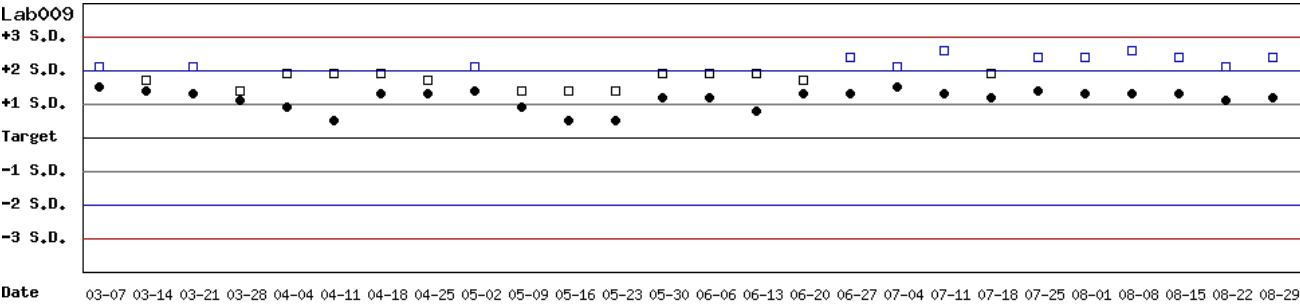


月統計 : 2025 08 Change ; 累計期間 (CUM) : 自 2025 03 07 至 2025 08 29 Change

TOP

Control N SDI 分佈圖

Lot No.: BA1116N ; 檢驗期間 : 2025-03-07 至 2025-08-29 ; 目標值 (Target) : 12.8 ; 參考標準差 (SD) : 1.28

單位 : Lab009

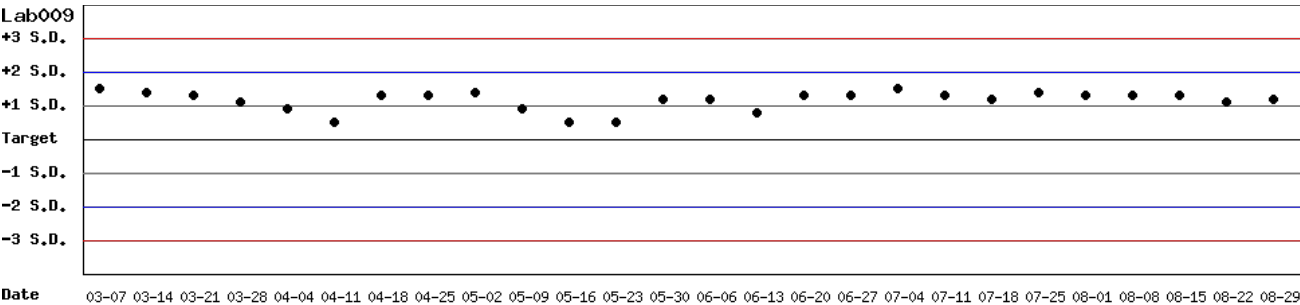


月統計 : 2025 08 Change ; 累計期間 (CUM) : 自 2025 03 07 至 2025 08 29 Change

TOP

Control D SDI 分佈圖

Lot No.: BA1116D; 檢驗期間: 2025-03-07 至 2025-08-29; 目標值 (Target): 4.2; 參考標準差 (SD): 0.42

單位: Lab009

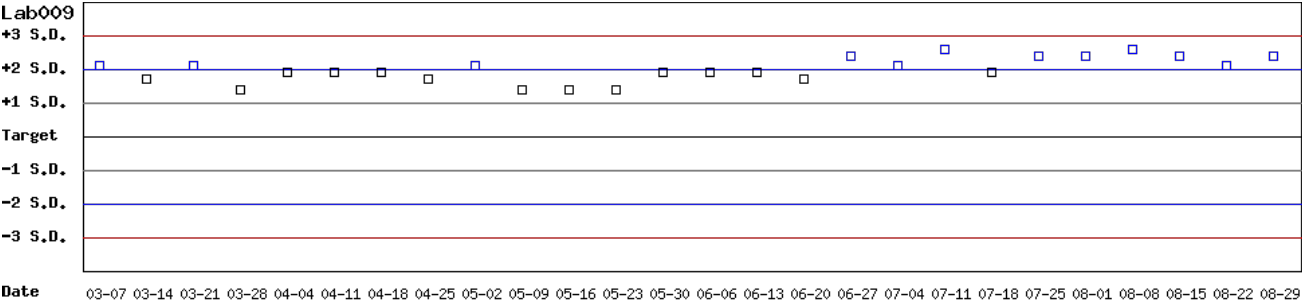

月統計: 2025 08 Change ; 累計期間 (CUM): 自 2025 03 07 至 2025 08 29 Change

[TOP](#)

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表一）

選擇 LotNo： BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit： 3 - Innovation Change

列印總表一

月統計表

統計月份： 2025 08 Change

		Control N (Lot No.: BA1116N)								Control D (Lot No.: BA1116D)							
UnitID 	Reagent Kit (Code) 	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab009	3	12.8	14.4	5	0.1	0.7	13.9	20	>6	4.2	5.2	5	0.1	1.9	27.7	20	-2.0
Lab051	3	14.4	18.2	4	0.3	1.6	29.7	20	-4.0	5.7	4.9	4	0.1	2.0	18.1	20	3.0
Total	-	-	16.1	9	2.0	12.4	-	-	-	-	5.0	9	0.2	4.0	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[\[TOP\]](#)

累計統計表

累計期間：自 2022 06 01 至 2025 08 31 Change

		Control N (Lot No.: BA1116N)								Control D (Lot No.: BA1116D)							
UnitID 	Reagent Kit (Code) 	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/ gHb)	Mean (U/ gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab004	3	14.8	14.8	143	0.4	2.7	5.4	30	>6	4.6	4.6	143	0.2	4.3	8.7	30	>6
Lab009	3	12.8	14.2	66	0.4	2.8	16.6	20	3.2	4.2	5.1	66	0.2	3.9	29.3	20	-0.4
Lab011	3	15.2	14.3	82	1.1	7.7	21.3	20	1.8	5.4	4.8	82	0.4	8.3	27.8	20	1.1
Lab013	3	15.1	15.3	176	0.3	2.0	5.2	20	>6	5.4	5.5	176	0.2	3.6	9.1	20	5.0
Lab014	3	14.4	14.3	244	0.3	2.1	4.9	20	>6	5.7	5.6	244	0.2	3.6	8.9	20	5.1
Lab016	3	12.8	12.8	92	0.4	3.1	6.3	20	>6	4.2	4.3	92	0.2	4.7	11.7	20	3.7
Lab029	3	14.8	15.1	123	0.6	4.0	10.0	20	4.5	5.1	5.3	123	0.3	5.7	15.2	20	2.8
Lab044	3	14.7	14.7	83	0.5	3.4	6.8	20	5.9	5.1	5.1	83	0.2	3.9	7.8	20	5.1
Lab051	3	14.4	15.7	92	2.1	13.4	35.8	20	0.8	5.7	5.2	92	1.4	26.9	62.6	20	0.4
Total	-	-	14.6	1101	1.1	7.5	-	-	-	-	5.1	1101	0.6	11.8	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Innovation	3

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表二）

選擇 LotNo：BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit：3 - Innovation Change

列印總表二

Control N 月統計與累計統計之比較表

Control N (Lot No.: BA1116N)																	
Month (2025/08)										CUM (2022/06/01~2025/08/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab004	3	14.8	-	0	-	-	-	30	-	14.8	14.8	143	0.4	2.7	5.4	30	>6
Lab009	3	12.8	14.4	5	0.1	0.7	13.9	20	>6	12.8	14.2	66	0.4	2.8	16.6	20	3.2
Lab011	3	15.2	-	0	-	-	-	20	-	15.2	14.3	82	1.1	7.7	21.3	20	1.8
Lab013	3	15.1	-	0	-	-	-	20	-	15.1	15.3	176	0.3	2.0	5.2	20	>6
Lab014	3	14.4	-	0	-	-	-	20	-	14.4	14.3	244	0.3	2.1	4.9	20	>6
Lab016	3	12.8	-	0	-	-	-	20	-	12.8	12.8	92	0.4	3.1	6.3	20	>6
Lab029	3	14.8	-	0	-	-	-	20	-	14.8	15.1	123	0.6	4.0	10.0	20	4.5
Lab044	3	14.7	-	0	-	-	-	20	-	14.7	14.7	83	0.5	3.4	6.8	20	5.9
Lab051	3	14.4	18.2	4	0.3	1.6	29.7	20	-4.0	14.4	15.7	92	2.1	13.4	35.8	20	0.8
Total	-	-	16.1	9	2.0	12.4	-	-	-	-	14.6	1101	1.1	7.5	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計：2025 08 Change 累計期間 (CUM)：自 2022 06 01 至 2025 08 31 Change

[\[TOP\]](#)

Control D 月統計與累計統計之比較表

Control D (Lot No.: BA1116D)																	
Month (2025/08)										CUM (2022/06/01~2025/08/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab004	3	4.6	-	0	-	-	-	30	-	4.6	4.6	143	0.2	4.3	8.7	30	>6
Lab009	3	4.2	5.2	5	0.1	1.9	27.7	20	-2.0	4.2	5.1	66	0.2	3.9	29.3	20	-0.4
Lab011	3	5.4	-	0	-	-	-	20	-	5.4	4.8	82	0.4	8.3	27.8	20	1.1
Lab013	3	5.4	-	0	-	-	-	20	-	5.4	5.5	176	0.2	3.6	9.1	20	5.0
Lab014	3	5.7	-	0	-	-	-	20	-	5.7	5.6	244	0.2	3.6	8.9	20	5.1
Lab016	3	4.2	-	0	-	-	-	20	-	4.2	4.3	92	0.2	4.7	11.7	20	3.7
Lab029	3	5.1	-	0	-	-	-	20	-	5.1	5.3	123	0.3	5.7	15.2	20	2.8
Lab044	3	5.1	-	0	-	-	-	20	-	5.1	5.1	83	0.2	3.9	7.8	20	5.1
Lab051	3	5.7	4.9	4	0.1	2.0	18.1	20	3.0	5.7	5.2	92	1.4	26.9	62.6	20	0.4
Total	-	-	5.0	9	0.2	4.0	-	-	-	-	5.1	1101	0.6	11.8	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計：2025 08 Change 累計期間 (CUM)：自 2022 06 01 至 2025 08 31 Change

[\[TOP\]](#)

Reagent Kit	Reagent Code
Innovation	3