


Summary Report of IQC program for G6PD Quantitative Test - AMP Group - June 2015 -

I. The statistic results of all laboratories in this month

G6PD	Control N (Lot No.:AC1203N)	Control D (Lot No.:AC1203D)
Labs	14	14
Received results number (n)	85	85
Median	12.6 (U/gHb)	3.1 (U/gHb)
Mean	12.7 (U/gHb)	3.1 (U/gHb)
SD	1.6	0.5
CV	12.6%	16.1%
Range of G6PD	9.4 ~ 15.6 (U/gHb)	2.2 ~ 4.5 (U/gHb)
Range of Hb	1.9 ~ 3.1 (g/dL)	1.5 ~ 5.1 (g/dL)

*The statistic results are calculated from all labs reported in this month

II. The distribution of G6PD reported for each lab in this survey

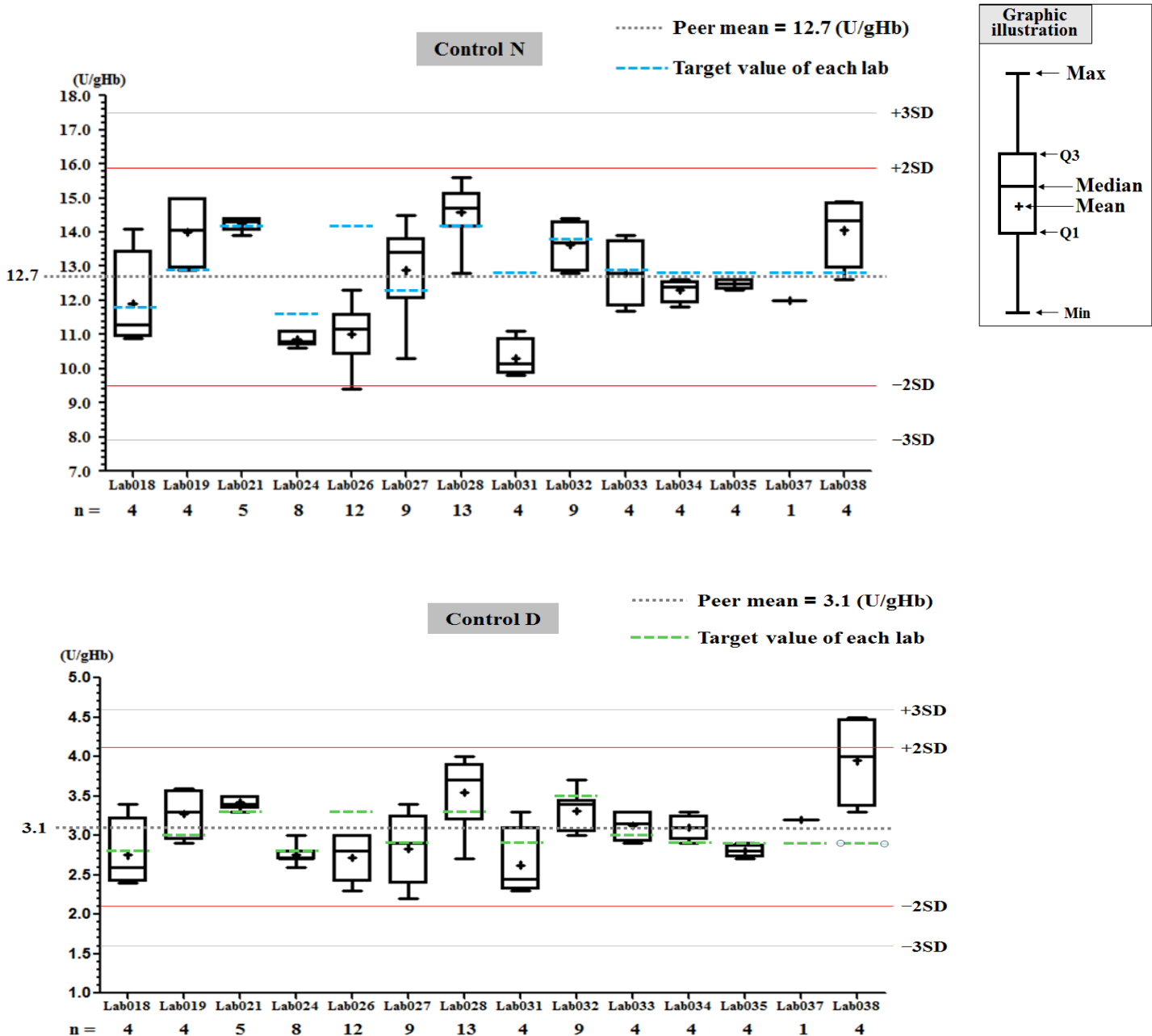

QC Chart of Internal Quality Control (IQC)for G6PD Quantitative Test

Select LotNo : AC1203N (2014-01-01 ~ 2100-12-31) [Change](#)

[Print Table](#)

Lab027

QC Control Lot No.	Control N		Control D	
	AC1203N		AC1203D	
Duration of the Analyzing	Month (2015/06)	CUM (2015/03/31~2015/06/30)	Month (2015/06)	CUM (2015/03/31~2015/06/30)
Runs (N)	9	26	9	26
Mean (U/gHb)	12.9	12.4	2.8	2.9
SD	1.3	1.8	0.4	0.5
CV (%)	10.1	14.5	14.3	17.2
Target Value (U/gHb)	12.3	12.3	2.9	2.9
Total Error (%)	25.0	29.8	32.0	34.5
TEa (%)	20	20	20	20
σ	1.5	1.3	1.2	1.2

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 x CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)


Month : 2015 06 [Change](#) ; Cumulative : from 2015 03 31 to 2015 06 30 [Change](#)

[TOP](#)

Control N SDI QC Chart

Lot No.: AC1203N ; Duration : 2015-03-31 to 2015-06-30 ; Target : 12.3 ; SD : 1.70

Lab027

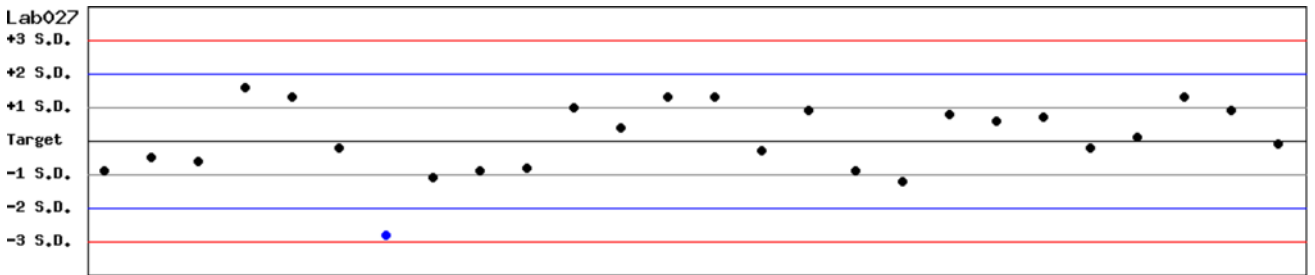

Date 03-31 04-07 04-09 04-14 04-16 04-21 04-23 04-28 04-30 05-05 05-07 05-12 05-14 05-19 05-21 05-26 05-28 06-02 06-04 06-09 06-11 06-16 06-18 06-23 06-25 06-30


Month : 2015 06 [Change](#) ; Cumulative : from 2015 03 31 to 2015 06 30 [Change](#)

[TOP](#)

Control D SDI QC Chart

Lot No.: AC1203D ; Duration : 2015-03-31 to 2015-06-30 ; Target : 2.9 ; SD : 0.50

Lab027

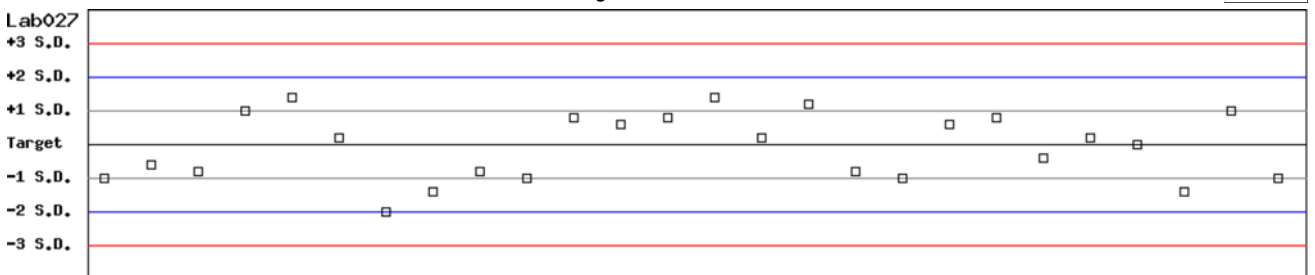

Date 03-31 04-07 04-09 04-14 04-16 04-21 04-23 04-28 04-30 05-05 05-07 05-12 05-14 05-19 05-21 05-26 05-28 06-02 06-04 06-09 06-11 06-16 06-18 06-23 06-25 06-30

Month : 2015 06 [Change](#) ; Cumulative : from 2015 03 31 to 2015 06 30 [Change](#)


[TOP](#)

Control N and Control D SDI QC Chart

Lot No.: AC1203N ; Duration : 2015-03-31 to 2015-06-30 ; Target : 12.3 ; SD : 1.70 (●)

Lot No.: AC1203D ; Duration : 2015-03-31 to 2015-06-30 ; Target : 2.9 ; SD : 0.50 (□)

Lab027

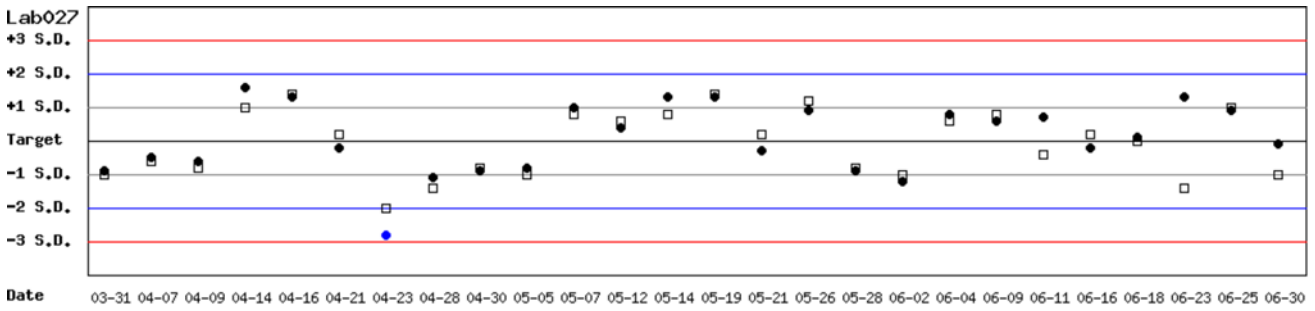

Month : 2015 06 ; Cumulative : from 2015 03 31 to 2015 06 30

[\[TOP\]](#)

Peer Group Statistics (Table 1)

Select LotNo : AC1203N (2014-01-01 ~ 2100-12-31) ▼ Change

Select Reagent Kit : 2 - AMP ▼ Change

[Print Table 1](#)

Monthly

Month : 2015 ▼ 06 ▼ Change

UnitID †	Reagent Kit (Code) †	Control N (Lot No.: AC1203N)								Control D (Lot No.: AC1203D)							
		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	2	11.8	11.9	4	1.5	12.6	26.1	20	1.5	2.8	2.8	4	0.5	17.9	35.7	20	1.1
Lab019	2	12.9	14.0	4	1.2	8.6	25.7	20	1.3	3.0	3.3	4	0.3	9.1	28.2	20	1.1
Lab021	2	14.2	14.2	5	0.2	1.4	2.8	20	>6	3.3	3.4	5	0.1	2.9	8.9	20	5.9
Lab024	2	11.6	10.9	8	0.2	1.8	9.7	20	>6	2.7	2.8	8	0.1	3.6	10.8	20	4.5
Lab026	2	14.2	11.0	12	0.8	7.3	37.1	20	-0.3	3.3	2.7	12	0.3	11.1	40.4	20	0.2
Lab027	2	12.3	12.9	9	1.3	10.1	25.0	20	1.5	2.9	2.8	9	0.4	14.3	32.0	20	1.2
Lab028	2	14.2	14.6	13	0.8	5.5	13.8	20	3.1	3.3	3.5	13	0.5	14.3	34.6	20	1.0
Lab031	2	12.8	10.3	4	0.6	5.8	31.2	20	0.1	2.9	2.6	4	0.5	19.2	48.8	20	0.5
Lab032	2	13.8	13.6	9	0.7	5.1	11.7	20	3.6	3.5	3.3	9	0.2	6.1	17.8	20	2.3
Lab033	2	12.9	12.8	4	1.0	7.8	16.4	20	2.5	3.0	3.1	4	0.2	6.5	16.2	20	2.6
Lab034	2	12.8	12.3	4	0.3	2.4	8.8	20	>6	2.9	3.1	4	0.2	6.5	19.8	20	2.0
Lab035	2	12.8	12.5	4	0.2	1.6	5.5	20	>6	2.9	2.8	4	0.1	3.6	10.6	20	4.6
Lab037	2	12.8	12.0	1	-	-	-	20	-	2.9	3.2	1	-	-	-	20	-
Lab038	2	12.8	14.1	4	1.0	7.1	24.3	20	1.4	2.9	4.0	4	0.6	15.0	67.9	20	-1.2
Total	-	-	12.7	85	1.6	12.6	-	-	-	-	3.1	85	0.5	16.1	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[TOP](#)

Cumulative

Cumulative : from 2014 ▼ 02 ▼ 01 ▼ to 2015 ▼ 06 ▼ 30 ▼ Change

UnitID †	Reagent Kit (Code) †	Control N (Lot No.: AC1203N)								Control D (Lot No.: AC1203D)							
		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab017	2	14.2	13.8	36	1.6	11.6	26.0	20	1.5	3.3	3.9	36	1.4	35.9	90.0	20	0.1
Lab018	2	11.8	10.5	49	1.7	16.2	43.4	20	0.6	2.8	2.4	49	1.2	50.0	114.3	20	0.1
Lab019	2	12.9	11.9	68	1.3	10.9	29.6	20	1.1	3.0	2.7	68	0.4	14.8	39.6	20	0.7
Lab021	2	14.2	14.2	14	0.2	1.4	2.8	20	>6	3.3	3.5	14	0.1	2.9	11.8	20	4.8
Lab024	2	11.6	11.1	39	1.3	11.7	27.7	20	1.3	2.7	2.7	39	0.4	14.8	29.6	20	1.4
Lab026	2	14.2	11.4	174	1.7	14.9	49.5	20	0.0	3.3	2.8	174	0.5	17.9	50.9	20	0.3
Lab027	2	12.3	11.7	104	1.7	14.5	33.9	20	1.0	2.9	2.7	104	0.5	18.5	43.9	20	0.7
Lab028	2	14.2	13.5	100	1.4	10.4	25.7	20	1.4	3.3	3.3	100	0.8	24.2	48.5	20	0.8
Lab031	2	12.8	9.9	45	1.4	14.1	50.9	20	-0.2	2.9	2.4	45	0.4	16.7	50.6	20	0.2
Lab032	2	13.8	10.7	133	1.7	15.9	54.2	20	-0.2	3.5	2.5	133	0.5	20.0	68.6	20	-0.4
Lab033	2	12.9	12.0	66	1.1	9.2	25.3	20	1.4	3.0	2.9	66	0.3	10.3	24.0	20	1.6
Lab034	2	12.8	12.6	18	0.4	3.2	7.9	20	5.8	2.9	3.1	18	0.2	6.5	19.8	20	2.0
Lab035	2	12.8	12.9	15	1.1	8.5	17.8	20	2.3	2.9	3.1	15	0.4	12.9	32.7	20	1.0
Lab037	2	12.8	12.9	3	0.9	7.0	14.7	20	2.7	2.9	3.2	3	0.2	6.3	22.8	20	1.5
Lab038	2	12.8	14.1	4	1.0	7.1	24.3	20	1.4	2.9	4.0	4	0.6	15.0	67.9	20	-1.2
Total	-	-	11.7	868	1.9	16.2	-	-	-	-	2.8	868	0.7	25.0	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[TOP](#)

Reagent Kit	Reagent Code
AMP	2

Peer Group Statistics (Table 2)

Select LotNo :

Select Reagent Kit :

Control N Month vs. Cumulative

		Control N (Lot No.: AC1203N)															
		Month (2015/06)								CUM (2014/02/01~2015/06/30)							
UnitID ↑	Reagent Kit (Code) ↑	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab017	2	14.2	-	0	-	-	-	20	-	14.2	13.8	36	1.6	11.6	26.0	20	1.5
Lab018	2	11.8	11.9	4	1.5	12.6	26.1	20	1.5	11.8	10.5	49	1.7	16.2	43.4	20	0.6
Lab019	2	12.9	14.0	4	1.2	8.6	25.7	20	1.3	12.9	11.9	68	1.3	10.9	29.6	20	1.1
Lab021	2	14.2	14.2	5	0.2	1.4	2.8	20	>6	14.2	14.2	14	0.2	1.4	2.8	20	>6
Lab024	2	11.6	10.9	8	0.2	1.8	9.7	20	>6	11.6	11.1	39	1.3	11.7	27.7	20	1.3
Lab026	2	14.2	11.0	12	0.8	7.3	37.1	20	-0.3	14.2	11.4	174	1.7	14.9	49.5	20	0.0
Lab027	2	12.3	12.9	9	1.3	10.1	25.0	20	1.5	12.3	11.7	104	1.7	14.5	33.9	20	1.0
Lab028	2	14.2	14.6	13	0.8	5.5	13.8	20	3.1	14.2	13.5	100	1.4	10.4	25.7	20	1.4
Lab031	2	12.8	10.3	4	0.6	5.8	31.2	20	0.1	12.8	9.9	45	1.4	14.1	50.9	20	-0.2
Lab032	2	13.8	13.6	9	0.7	5.1	11.7	20	3.6	13.8	10.7	133	1.7	15.9	54.2	20	-0.2
Lab033	2	12.9	12.8	4	1.0	7.8	16.4	20	2.5	12.9	12.0	66	1.1	9.2	25.3	20	1.4
Lab034	2	12.8	12.3	4	0.3	2.4	8.8	20	>6	12.8	12.6	18	0.4	3.2	7.9	20	5.8
Lab035	2	12.8	12.5	4	0.2	1.6	5.5	20	>6	12.8	12.9	15	1.1	8.5	17.8	20	2.3
Lab037	2	12.8	12.0	1	-	-	-	20	-	12.8	12.9	3	0.9	7.0	14.7	20	2.7
Lab038	2	12.8	14.1	4	1.0	7.1	24.3	20	1.4	12.8	14.1	4	1.0	7.1	24.3	20	1.4
Total	-	-	12.7	85	1.6	12.6	-	-	-	-	11.7	868	1.9	16.2	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%
 TE : Total Error(%) = Bias (%) + 2 x CV (%)
 σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month :

Cumulative : from to

[TOP](#)

Control D Month vs. Cumulative

		Control D (Lot No.: AC1203D)															
		Month (2015/06)								CUM (2014/02/01~2015/06/30)							
UnitID ↑	Reagent Kit (Code) ↑	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab017	2	3.3	-	0	-	-	-	20	-	3.3	3.9	36	1.4	35.9	90.0	20	0.1
Lab018	2	2.8	2.8	4	0.5	17.9	35.7	20	1.1	2.8	2.4	49	1.2	50.0	114.3	20	0.1
Lab019	2	3.0	3.3	4	0.3	9.1	28.2	20	1.1	3.0	2.7	68	0.4	14.8	39.6	20	0.7
Lab021	2	3.3	3.4	5	0.1	2.9	8.9	20	5.9	3.3	3.5	14	0.1	2.9	11.8	20	4.8
Lab024	2	2.7	2.8	8	0.1	3.6	10.8	20	4.5	2.7	2.7	39	0.4	14.8	29.6	20	1.4
Lab026	2	3.3	2.7	12	0.3	11.1	40.4	20	0.2	3.3	2.8	174	0.5	17.9	50.9	20	0.3
Lab027	2	2.9	2.8	9	0.4	14.3	32.0	20	1.2	2.9	2.7	104	0.5	18.5	43.9	20	0.7
Lab028	2	3.3	3.5	13	0.5	14.3	34.6	20	1.0	3.3	3.3	100	0.8	24.2	48.5	20	0.8
Lab031	2	2.9	2.6	4	0.5	19.2	48.8	20	0.5	2.9	2.4	45	0.4	16.7	50.6	20	0.2
Lab032	2	3.5	3.3	9	0.2	6.1	17.8	20	2.3	3.5	2.5	133	0.5	20.0	68.6	20	-0.4
Lab033	2	3.0	3.1	4	0.2	6.5	16.2	20	2.6	3.0	2.9	66	0.3	10.3	24.0	20	1.6
Lab034	2	2.9	3.1	4	0.2	6.5	19.8	20	2.0	2.9	3.1	18	0.2	6.5	19.8	20	2.0
Lab035	2	2.9	2.8	4	0.1	3.6	10.6	20	4.6	2.9	3.1	15	0.4	12.9	32.7	20	1.0
Lab037	2	2.9	3.2	1	-	-	-	20	-	2.9	3.2	3	0.2	6.3	22.8	20	1.5
Lab038	2	2.9	4.0	4	0.6	15.0	67.9	20	-1.2	2.9	4.0	4	0.6	15.0	67.9	20	-1.2
Total	-	-	3.1	85	0.5	16.1	-	-	-	-	2.8	868	0.7	25.0	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%
 TE : Total Error(%) = Bias (%) + 2 x CV (%)
 σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month :

Cumulative : from to

[TOP](#)

Reagent Kit	Reagent Code
AMP	2