

G6PD 定量檢驗內部品管統計結果及 SDI 分佈圖

選擇 LotNo : AH1220N (2019-01-28 ~ 2100-12-31) Change

列印圖表

檢驗單位 : Lab044

| QC Control Lot No. | Control N | | Control D | |
|--------------------|-----------------|-----------------------------|-----------------|-----------------------------|
| | AH1220N | | AH1220D | |
| 檢驗期間 | Month (2022/08) | CUM (2022/03/04~2022/08/31) | Month (2022/08) | CUM (2022/03/04~2022/08/31) |
| 累計點數 (N) | 4 | 26 | 4 | 26 |
| 累計平均值 (Mean) U/gHb | 16.9 | 17.4 | 5.5 | 5.5 |
| 累計標準差 (SD) | 0.1 | 0.4 | 0.1 | 0.2 |
| 累計變異係數 (CV) % | 0.6 | 2.3 | 1.8 | 3.6 |
| 目標值 (Target) U/gHb | 17.7 | 17.7 | 5.6 | 5.6 |
| 總誤差 (TE) % | 5.7 | 6.3 | 5.4 | 9.1 |
| 可容許總誤差 (TEa) % | 20 | 20 | 20 | 20 |
| σ | >6 | >6 | >6 | 5.1 |

誤差 (Bias)% = $\{ [| \text{累計平均值 (Mean)} - \text{目標值 (Target)} |] / \text{目標值 (Target)} \} \times 100\%$
 總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)
 σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2022 08 Change ; 累計期間 (CUM) : 自 2022 03 04 至 2022 08 31 Change


[TOP](#)

Control N 及 Control D SDI 分佈圖

Lot No.: AH1220N; 檢驗期間: 2022-03-04 至 2022-08-31; 目標值 (Target): 17.7; 參考標準差 (SD): 0.90 (●)

Lot No.: AH1220D; 檢驗期間: 2022-03-04 至 2022-08-31; 目標值 (Target): 5.6; 參考標準差 (SD): 0.30 (□)

單位: Lab044

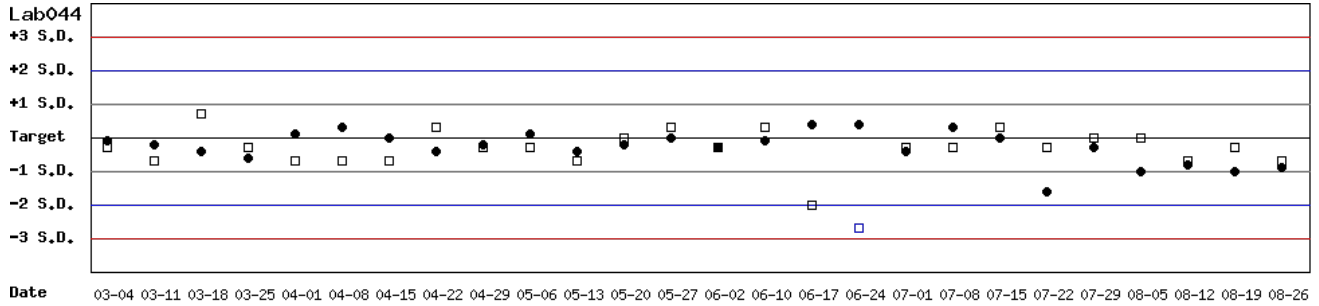


月統計 : 2022 08 Change ; 累計期間 (CUM) : 自 2022 03 04 至 2022 08 31 Change

[TOP](#)

Control N SDI 分佈圖

Lot No.: AH1220N; 檢驗期間: 2022-03-04 至 2022-08-31; 目標值 (Target): 17.7; 參考標準差 (SD): 0.90

單位: Lab044

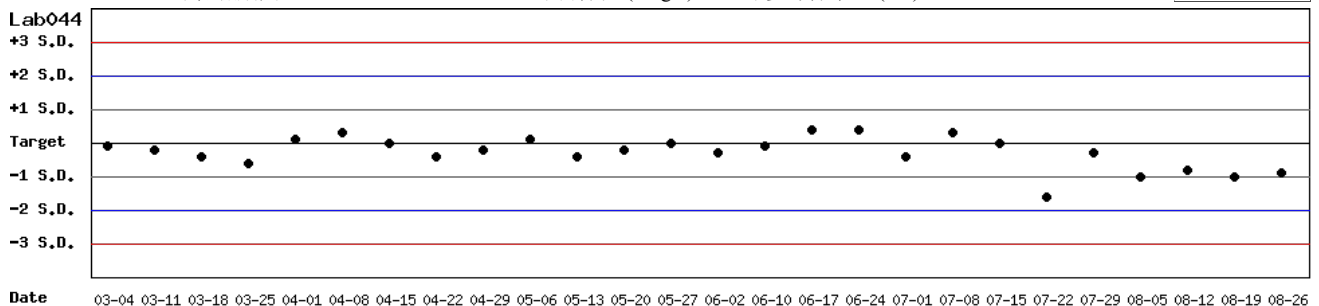


月統計 : 2022 08 Change ; 累計期間 (CUM) : 自 2022 03 04 至 2022 08 31 Change

[TOP](#)

Control D SDI 分佈圖

Lot No.: AH1220D ; 檢驗期間 : 2022-03-04 至 2022-08-31 ; 目標值 (Target) : 5.6 ; 參考標準差 (SD) : 0.30

單位 : Lab044

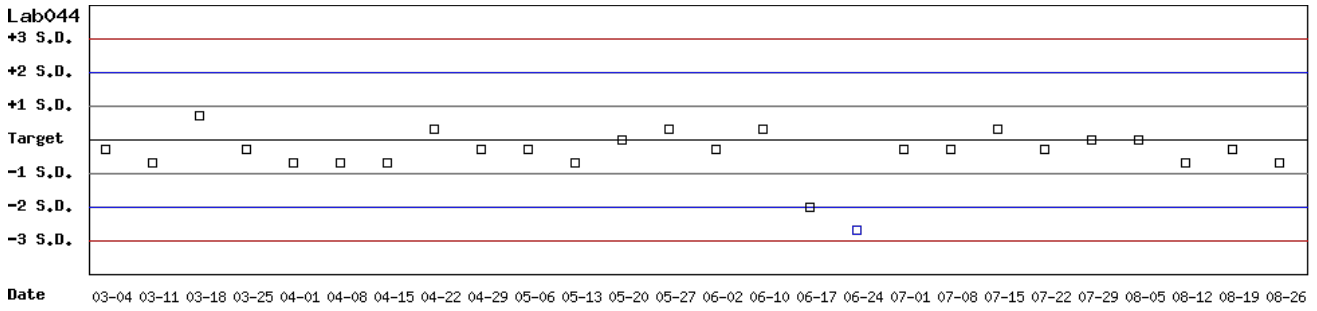

月統計 : ; 累計期間 (CUM) : 自 至

[TOP](#)

各檢驗單位 G6PD 定量檢驗 內部品管統計結果 (總表一)

選擇 LotNo : AH1220N (2019-01-28 ~ 2100-12-31) Change

選擇 Reagent Kit : 3 - Innovation Change

[列印總表一](#)

月統計表

統計月份 : 2022 08 Change

| UnitID ↑ | Reagent Kit (Code) ↑ | Control N (Lot No.: AH1220N) | | | | | | | | Control D (Lot No.: AH1220D) | | | | | | | |
|------------------------|-----------------------------------|------------------------------|--------------|------------|-----|--------|--------|---------|----|------------------------------|--------------|------------|-----|--------|--------|---------|----|
| | | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab044 | 3 | 17.7 | 16.9 | 4 | 0.1 | 0.6 | 5.7 | 20 | >6 | 5.6 | 5.5 | 4 | 0.1 | 1.8 | 5.4 | 20 | >6 |
| Total | - | - | 16.9 | 4 | 0.1 | 0.6 | - | - | - | - | 5.5 | 4 | 0.1 | 1.8 | - | - | - |

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[TOP](#)

累計統計表

累計期間 : 自 2019 02 01 至 2022 08 31 Change

| UnitID ↑ | Reagent Kit (Code) ↑ | Control N (Lot No.: AH1220N) | | | | | | | | Control D (Lot No.: AH1220D) | | | | | | | |
|------------------------|-----------------------------------|------------------------------|--------------|------------|-----|--------|--------|---------|-----|------------------------------|--------------|------------|-----|--------|--------|---------|-----|
| | | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab004 | 3 | 17.1 | 17.3 | 99 | 0.4 | 2.3 | 5.8 | 20 | >6 | 5.6 | 5.5 | 99 | 0.2 | 3.6 | 9.1 | 20 | 5.1 |
| Lab013 | 3 | 15.8 | 15.9 | 213 | 0.2 | 1.3 | 3.1 | 20 | >6 | 5.2 | 5.1 | 213 | 0.6 | 11.8 | 25.5 | 20 | 1.5 |
| Lab014 | 3 | 16.6 | 16.2 | 124 | 0.7 | 4.3 | 11.1 | 20 | 4.1 | 5.4 | 5.4 | 124 | 0.2 | 3.7 | 7.4 | 20 | 5.4 |
| Lab016 | 3 | 15.3 | 16.1 | 78 | 0.2 | 1.2 | 7.7 | 20 | >6 | 5.2 | 5.2 | 78 | 0.1 | 1.9 | 3.8 | 20 | >6 |
| Lab029 | 3 | 16.2 | 16.2 | 60 | 0.5 | 3.1 | 6.2 | 20 | >6 | 5.3 | 5.4 | 60 | 0.3 | 5.6 | 13.0 | 20 | 3.2 |
| Lab044 | 3 | 17.7 | 17.1 | 251 | 0.7 | 4.1 | 11.6 | 20 | 4.1 | 5.6 | 5.4 | 251 | 0.2 | 3.7 | 11.0 | 20 | 4.4 |
| Total | - | - | 16.5 | 825 | 0.7 | 4.2 | - | - | - | - | 5.3 | 825 | 0.4 | 7.5 | - | - | - |

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

[TOP](#)

| Reagent Kit | Reagent Code |
|-------------|--------------|
| Innovation | 3 |

各檢驗單位 G6PD 定量檢驗 內部品管統計結果 (總表二)

選擇 LotNo : AH1220N (2019-01-28 ~ 2100-12-31) Change

選擇 Reagent Kit : 3 - Innovation Change

列印總表二

Control N 月統計與累計統計之比較表

| | | Control N (Lot No.: AH1220N) | | | | | | | | | | | | | | | |
|------------------------|-----------------------------------|------------------------------|--------------|------------|-----|--------|--------|---------|----|-----------------------------|--------------|------------|-----|--------|--------|---------|-----|
| | | Month (2022/08) | | | | | | | | CUM (2019/02/01~2022/08/31) | | | | | | | |
| UnitID ↑ | Reagent Kit (Code) ↑ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab004 | 3 | 17.1 | - | 0 | - | - | - | 20 | - | 17.1 | 17.3 | 99 | 0.4 | 2.3 | 5.8 | 20 | >6 |
| Lab013 | 3 | 15.8 | - | 0 | - | - | - | 20 | - | 15.8 | 15.9 | 213 | 0.2 | 1.3 | 3.1 | 20 | >6 |
| Lab014 | 3 | 16.6 | - | 0 | - | - | - | 20 | - | 16.6 | 16.2 | 124 | 0.7 | 4.3 | 11.1 | 20 | 4.1 |
| Lab016 | 3 | 15.3 | - | 0 | - | - | - | 20 | - | 15.3 | 16.1 | 78 | 0.2 | 1.2 | 7.7 | 20 | >6 |
| Lab029 | 3 | 16.2 | - | 0 | - | - | - | 20 | - | 16.2 | 16.2 | 60 | 0.5 | 3.1 | 6.2 | 20 | >6 |
| Lab044 | 3 | 17.7 | 16.9 | 4 | 0.1 | 0.6 | 5.7 | 20 | >6 | 17.7 | 17.1 | 251 | 0.7 | 4.1 | 11.6 | 20 | 4.1 |
| Total | - | - | 16.9 | 4 | 0.1 | 0.6 | - | - | - | - | 16.5 | 825 | 0.7 | 4.2 | - | - | - |

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2022 ▼ 08 ▼ Change 累計期間 (CUM) : 自 2019 ▼ 02 ▼ 01 ▼ 至 2022 ▼ 08 ▼ 31 ▼ Change

[TOP](#)

Control D 月統計與累計統計之比較表

| | | Control D (Lot No.: AH1220D) | | | | | | | | | | | | | | | |
|------------------------|-----------------------------------|------------------------------|--------------|------------|-----|--------|--------|---------|----|-----------------------------|--------------|------------|-----|--------|--------|---------|-----|
| | | Month (2022/08) | | | | | | | | CUM (2019/02/01~2022/08/31) | | | | | | | |
| UnitID ↑ | Reagent Kit (Code) ↑ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ | Target (U/gHb) | Mean (U/gHb) | n for Mean | SD | CV (%) | TE (%) | TEa (%) | σ |
| Lab004 | 3 | 5.6 | - | 0 | - | - | - | 20 | - | 5.6 | 5.5 | 99 | 0.2 | 3.6 | 9.1 | 20 | 5.1 |
| Lab013 | 3 | 5.2 | - | 0 | - | - | - | 20 | - | 5.2 | 5.1 | 213 | 0.6 | 11.8 | 25.5 | 20 | 1.5 |
| Lab014 | 3 | 5.4 | - | 0 | - | - | - | 20 | - | 5.4 | 5.4 | 124 | 0.2 | 3.7 | 7.4 | 20 | 5.4 |
| Lab016 | 3 | 5.2 | - | 0 | - | - | - | 20 | - | 5.2 | 5.2 | 78 | 0.1 | 1.9 | 3.8 | 20 | >6 |
| Lab029 | 3 | 5.3 | - | 0 | - | - | - | 20 | - | 5.3 | 5.4 | 60 | 0.3 | 5.6 | 13.0 | 20 | 3.2 |
| Lab044 | 3 | 5.6 | 5.5 | 4 | 0.1 | 1.8 | 5.4 | 20 | >6 | 5.6 | 5.4 | 251 | 0.2 | 3.7 | 11.0 | 20 | 4.4 |
| Total | - | - | 5.5 | 4 | 0.1 | 1.8 | - | - | - | - | 5.3 | 825 | 0.4 | 7.5 | - | - | - |

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2022 ▼ 08 ▼ Change 累計期間 (CUM) : 自 2019 ▼ 02 ▼ 01 ▼ 至 2022 ▼ 08 ▼ 31 ▼ Change

[TOP](#)

| Reagent Kit | Reagent Code |
|-------------|--------------|
| Innovation | 3 |