

Summary Report of IQC program for G6PD Quantitative Test - Medicon Group - March 2025 -


I. The statistic results of all laboratories in this month

G6PD	Control N (Lot No.:BJ0922N)	Control D (Lot No.: BJ0922D)
Labs	6	6
Received results number (n)	24	24
Median	14.7 (U/g Hb)	4.8 (U/g Hb)
Mean	14.6 (U/g Hb)	4.9 (U/g Hb)
SD	0.5	0.5
CV	3.4%	10.2%
Range of G6PD	13.2 ~ 15.6 (U/g Hb)	4.0 ~ 5.7 (U/g Hb)
Range of Hb	2.0 ~ 2.3 (g/dL)	1.8 ~ 2.2 (g/dL)

*The statistic results are calculated from all labs reported in this month

** G6PD Method = Medicon reagent kit, 37°C


II. The distribution of G6PD reported for each lab in this survey

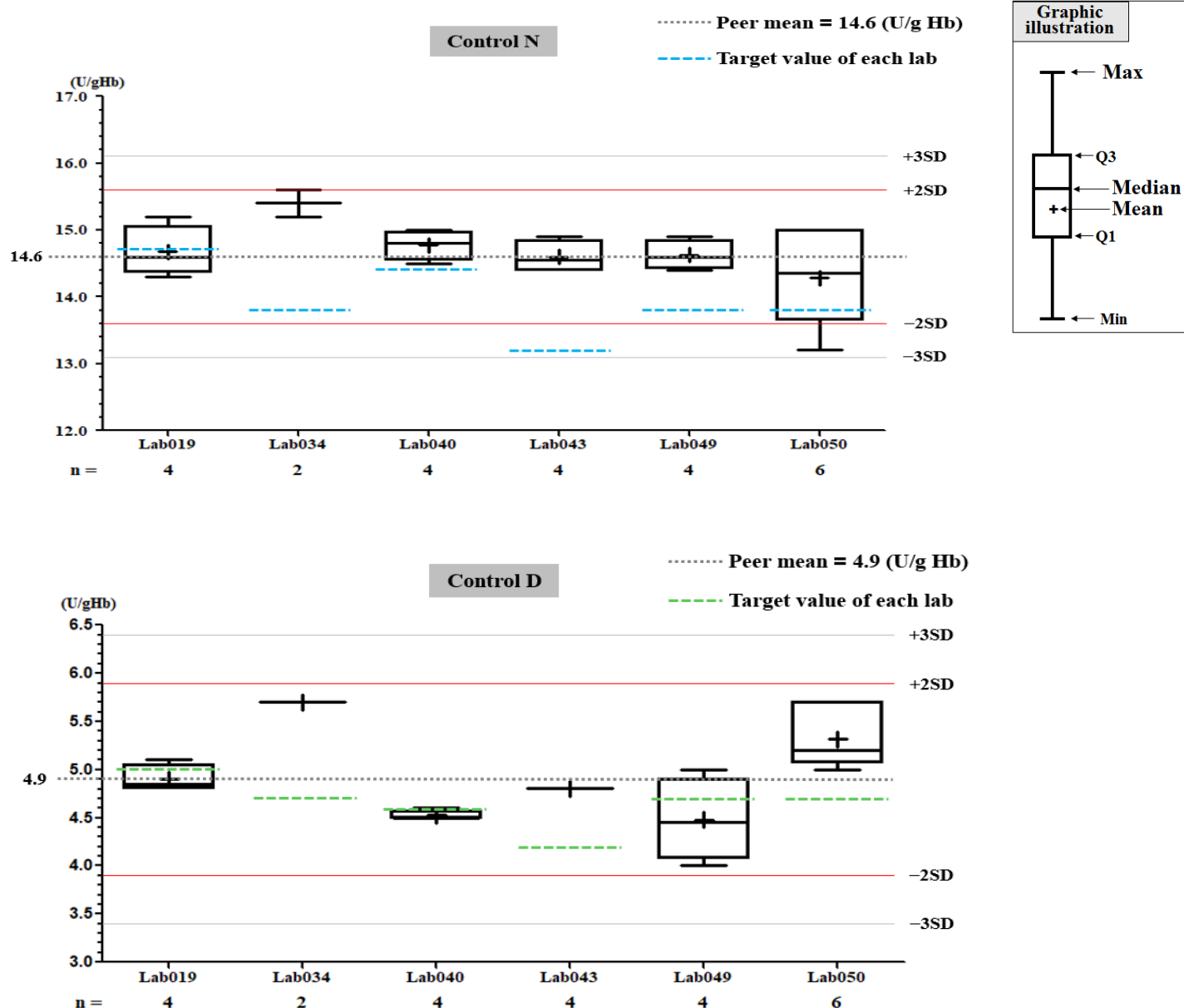

Control D SDI QC Chart

Lot No.: BJ0922D ; Duration : 2024-10-01 to 2025-03-31 ; Target : 5.0 ; SD : 0.30

Lab019

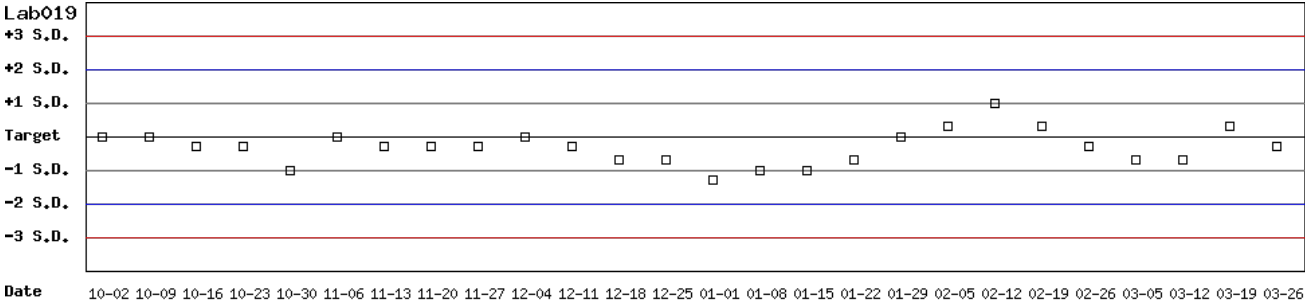

Month : 2025 03 Change ; Cumulative : from 2024 10 01 to 2025 03 31 Change

[\[TOP\]](#)

Peer Group Statistics (Table 1)

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

Print Table 1

Monthly

Month : 2025 03 Change

		Control N (Lot No.: BJ0922N)									Control D (Lot No.: BJ0922D)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab019	5	14.7	14.7	4	0.4	2.7	5.4	20	>6		5.0	4.9	4	0.1	2.0	6.1	20	>6
Lab034	5	13.8	15.4	2	-	-	-	20	-		4.7	5.7	2	-	-	-	20	-
Lab040	5	14.4	14.8	4	0.2	1.4	5.5	20	>6		4.6	4.5	4	0.1	2.2	6.6	20	>6
Lab043	5	13.2	14.6	4	0.2	1.4	13.3	20	>6		4.2	4.8	4	0.0	0.0	14.3	20	>6
Lab049	5	13.8	14.6	4	0.2	1.4	8.5	20	>6		4.7	4.5	4	0.4	8.9	22.0	20	1.8
Lab050	5	13.8	14.3	6	0.7	4.9	13.4	20	3.3		4.7	5.3	6	0.3	5.7	24.1	20	1.3
Total	-	-	14.6	24	0.5	3.4	-	-	-		-	4.9	24	0.5	10.2	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 x CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Cumulative

Cumulative : from 2021 01 01 to 2025 03 31 Change

		Control N (Lot No.: BJ0922N)									Control D (Lot No.: BJ0922D)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	16.2	15.3	113	0.4	2.6	10.8	20	5.6		5.5	5.4	113	0.1	1.9	5.5	20	>6
Lab019	5	14.7	14.6	157	0.5	3.4	7.5	20	5.7		5.0	4.8	157	0.2	4.2	12.3	20	3.8
Lab020	5	14.4	14.4	90	0.5	3.5	6.9	20	5.7		4.6	4.7	90	0.2	4.3	10.7	20	4.1
Lab022	5	14.4	14.4	27	0.5	3.5	6.9	20	5.7		4.6	4.7	27	0.3	6.4	14.9	20	2.8
Lab025	5	14.4	14.2	13	0.5	3.5	8.4	20	5.3		5.4	5.1	13	0.3	5.9	17.3	20	2.4
Lab028	5	15.2	15.0	66	0.9	6.0	13.3	20	3.1		4.9	5.1	66	0.3	5.9	15.8	20	2.7
Lab032	5	14.9	14.3	153	0.4	2.8	9.6	20	5.7		4.8	4.5	153	0.2	4.4	15.1	20	3.1
Lab033	5	13.8	13.6	81	0.7	5.1	11.7	20	3.6		4.7	4.6	81	0.3	6.5	15.2	20	2.7
Lab034	5	13.8	14.4	33	0.5	3.5	11.3	20	4.5		4.7	5.2	33	0.2	3.8	18.3	20	2.5
Lab037	5	13.8	14.5	59	0.5	3.4	12.0	20	4.4		4.7	4.7	59	0.2	4.3	8.5	20	4.7
Lab040	5	14.4	14.8	152	0.2	1.4	5.5	20	>6		4.6	4.6	152	0.1	2.2	4.3	20	>6
Lab043	5	13.2	14.0	183	1.0	7.1	20.3	20	2.0		4.2	4.6	183	0.4	8.7	26.9	20	1.2
Lab049	5	13.8	14.4	251	0.6	4.2	12.7	20	3.7		4.7	4.7	251	0.3	6.4	12.8	20	3.1
Lab050	5	13.8	14.5	80	0.6	4.1	13.3	20	3.6		4.7	5.0	80	0.5	10.0	26.4	20	1.4
Total	-	-	14.5	1458	0.7	4.8	-	-	-		-	4.8	1458	0.4	8.3	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 x CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Medicon	5

Peer Group Statistics (Table 2)

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

[Print Table 2](#)

Control N Month vs. Cumulative

Control N (Lot No.: BJ0922N)																	
		Month (2025/03)								CUM (2021/01/01~2025/03/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	16.2	-	0	-	-	-	20	-	16.2	15.3	113	0.4	2.6	10.8	20	5.6
Lab019	5	14.7	14.7	4	0.4	2.7	5.4	20	>6	14.7	14.6	157	0.5	3.4	7.5	20	5.7
Lab020	5	14.4	-	0	-	-	-	20	-	14.4	14.4	90	0.5	3.5	6.9	20	5.7
Lab022	5	14.4	-	0	-	-	-	20	-	14.4	14.4	27	0.5	3.5	6.9	20	5.7
Lab025	5	14.4	-	0	-	-	-	20	-	14.4	14.2	13	0.5	3.5	8.4	20	5.3
Lab028	5	15.2	-	0	-	-	-	20	-	15.2	15.0	66	0.9	6.0	13.3	20	3.1
Lab032	5	14.9	-	0	-	-	-	20	-	14.9	14.3	153	0.4	2.8	9.6	20	5.7
Lab033	5	13.8	-	0	-	-	-	20	-	13.8	13.6	81	0.7	5.1	11.7	20	3.6
Lab034	5	13.8	15.4	2	-	-	-	20	-	13.8	14.4	33	0.5	3.5	11.3	20	4.5
Lab037	5	13.8	-	0	-	-	-	20	-	13.8	14.5	59	0.5	3.4	12.0	20	4.4
Lab040	5	14.4	14.8	4	0.2	1.4	5.5	20	>6	14.4	14.8	152	0.2	1.4	5.5	20	>6
Lab043	5	13.2	14.6	4	0.2	1.4	13.3	20	>6	13.2	14.0	183	1.0	7.1	20.3	20	2.0
Lab049	5	13.8	14.6	4	0.2	1.4	8.5	20	>6	13.8	14.4	251	0.6	4.2	12.7	20	3.7
Lab050	5	13.8	14.3	6	0.7	4.9	13.4	20	3.3	13.8	14.5	80	0.6	4.1	13.3	20	3.6
Total	-	-	14.6	24	0.5	3.4	-	-	-	-	14.5	1458	0.7	4.8	-	-	-

Bias (%) = $\left[\left(\left| \text{Mean} - \text{Target} \right| \right) / \text{Target} \right] \times 100\%$

TE : Total Error(%) = Bias (%) + $2 \times \text{CV}$ (%)

σ (Sigma) = $[\text{TEa}\% - \text{Bias}\%] / \text{CV}\%$

Month : 2025  03  Change

Cumulative : from 2021  01  01  to 2025  03  31  Change

[\[TOP\]](#)

Control D Month vs. Cumulative

Control D (Lot No.: BJ0922D)																	
		Month (2025/03)								CUM (2021/01/01~2025/03/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	5.5	-	0	-	-	-	20	-	5.5	5.4	113	0.1	1.9	5.5	20	>6
Lab019	5	5.0	4.9	4	0.1	2.0	6.1	20	>6	5.0	4.8	157	0.2	4.2	12.3	20	3.8
Lab020	5	4.6	-	0	-	-	-	20	-	4.6	4.7	90	0.2	4.3	10.7	20	4.1
Lab022	5	4.6	-	0	-	-	-	20	-	4.6	4.7	27	0.3	6.4	14.9	20	2.8
Lab025	5	5.4	-	0	-	-	-	20	-	5.4	5.1	13	0.3	5.9	17.3	20	2.4
Lab028	5	4.9	-	0	-	-	-	20	-	4.9	5.1	66	0.3	5.9	15.8	20	2.7
Lab032	5	4.8	-	0	-	-	-	20	-	4.8	4.5	153	0.2	4.4	15.1	20	3.1
Lab033	5	4.7	-	0	-	-	-	20	-	4.7	4.6	81	0.3	6.5	15.2	20	2.7
Lab034	5	4.7	5.7	2	-	-	-	20	-	4.7	5.2	33	0.2	3.8	18.3	20	2.5
Lab037	5	4.7	-	0	-	-	-	20	-	4.7	4.7	59	0.2	4.3	8.5	20	4.7
Lab040	5	4.6	4.5	4	0.1	2.2	6.6	20	>6	4.6	4.6	152	0.1	2.2	4.3	20	>6
Lab043	5	4.2	4.8	4	0.0	0.0	14.3	20	>6	4.2	4.6	183	0.4	8.7	26.9	20	1.2
Lab049	5	4.7	4.5	4	0.4	8.9	22.0	20	1.8	4.7	4.7	251	0.3	6.4	12.8	20	3.1
Lab050	5	4.7	5.3	6	0.3	5.7	24.1	20	1.3	4.7	5.0	80	0.5	10.0	26.4	20	1.4
Total	-	-	4.9	24	0.5	10.2	-	-	-	-	4.8	1458	0.4	8.3	-	-	-

Bias (%) = $\left[\left(\left| \text{Mean} - \text{Target} \right| \right) / \text{Target} \right] \times 100\%$

TE : Total Error(%) = Bias (%) + $2 \times \text{CV}$ (%)

σ (Sigma) = $[\text{TEa}\% - \text{Bias}\%] / \text{CV}\%$

Month : 2025  03  Change

Cumulative : from 2021  01  01  to 2025  03  31  Change

[\[TOP\]](#)

Reagent Kit	Reagent Code
Medicon	5