

G6PD 定量檢驗內部品管統計結果及 SDI 分佈圖

選擇 LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

列印圖表

檢驗單位 : Lab051

QC Control Lot No.	Control N		Control D	
	BA1116N		BA1116D	
檢驗期間	Month (2025/06)	CUM (2025/01/07~2025/06/27)	Month (2025/06)	CUM (2025/01/07~2025/06/27)
累計點數 (N)	3	26	3	26
累計平均值 (Mean) U/gHb	17.8	17.3	4.8	5.4
累計標準差 (SD)	0.3	2.7	0.2	2.6
累計變異係數 (CV) %	1.7	15.6	4.2	48.1
目標值 (Target) U/gHb	14.4	14.4	5.7	5.7
總誤差 (TE) %	27.0	51.4	24.1	101.6
可容許總誤差 (TEa) %	20	20	20	20
σ	-2.1	0.0	1.0	0.3

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 × 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計 : 2025 06 Change ; 累計期間 (CUM) : 自 2025 01 07 至 2025 06 27 Change


TOP

Control N 及 Control D SDI
分佈圖

Lot No.: BA1116N ; 檢驗期間 : 2025-01-07 至 2025-06-27 ; 目標值 (Target) : 14.4 ; 參考標準差 (SD) : 0.50 (●)

Lot No.: BA1116D ; 檢驗期間 : 2025-01-07 至 2025-06-27 ; 目標值 (Target) : 5.7 ; 參考標準差 (SD) : 0.20 (□)

單位 : Lab051

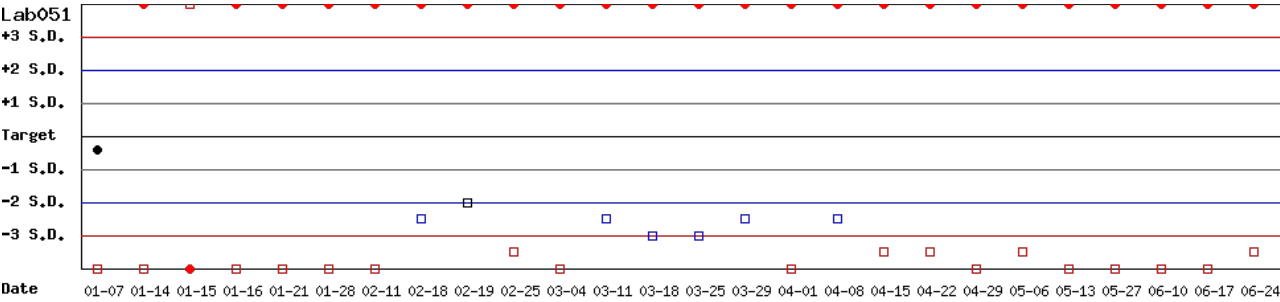


月統計 : 2025 06 Change ; 累計期間 (CUM) : 自 2025 01 07 至 2025 06 27 Change

TOP

Control N SDI 分佈圖

Lot No.: BA1116N ; 檢驗期間 : 2025-01-07 至 2025-06-27 ; 目標值 (Target) : 14.4 ; 參考標準差 (SD) : 0.50

單位 : Lab051

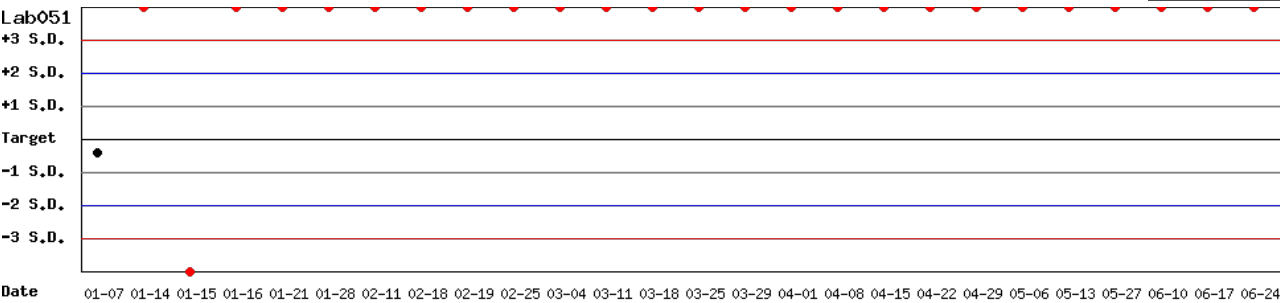


月統計 : 2025 06 Change ; 累計期間 (CUM) : 自 2025 01 07 至 2025 06 27 Change

TOP

Control D SDI 分佈圖

Lot No.: BA1116D ; 檢驗期間 : 2025-01-07 至 2025-06-27 ; 目標值 (Target) : 5.7 ; 參考標準差 (SD) : 0.20

單位 : Lab051

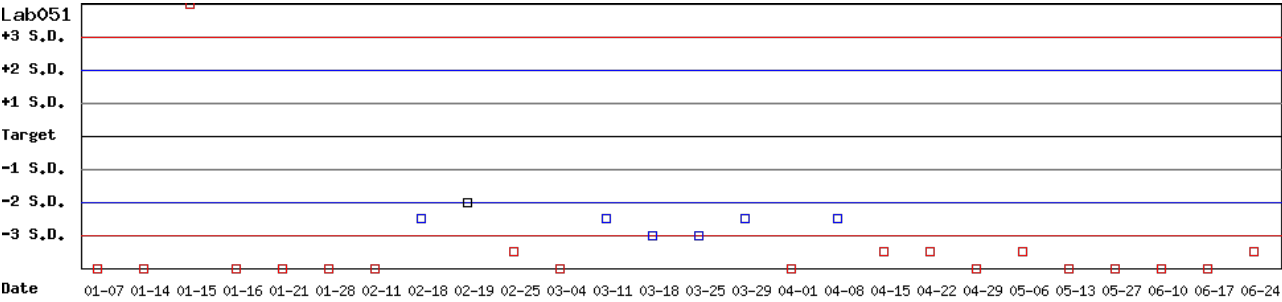

月統計 : 2025 ▾ 06 ▾ Change ; 累計期間 (CUM) : 自 2025 ▾ 01 ▾ 07 ▾ 至 2025 ▾ 06 ▾ 27 ▾ Change

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表一）

選擇 LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit : 3 - Innovation Change

列印總表一

月統計表

統計月份 : 2025 06 Change

		Control N (Lot No.: BA1116N)									Control D (Lot No.: BA1116D)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab009	3	12.8	14.3	4	0.3	2.1	15.9	20	3.9		4.2	5.0	4	0.1	2.0	23.0	20	0.5
Lab044	3	14.7	15.2	4	0.1	0.7	4.7	20	>6		5.1	5.2	4	0.1	1.9	5.8	20	>6
Lab051	3	14.4	17.8	3	0.3	1.7	27.0	20	-2.1		5.7	4.8	3	0.2	4.2	24.1	20	1.0
Total	-	-	15.6	11	1.5	9.6	-	-	-		-	5.0	11	0.2	4.0	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

TOP

累計統計表

累計期間 : 自 2022 06 01 至 2025 06 30 Change

		Control N (Lot No.: BA1116N)									Control D (Lot No.: BA1116D)								
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	
Lab004	3	14.8	14.8	143	0.4	2.7	5.4	30	>6		4.6	4.6	143	0.2	4.3	8.7	30	>6	
Lab009	3	12.8	14.1	57	0.4	2.8	15.8	20	3.5		4.2	5.1	57	0.2	3.9	29.3	20	-0.4	
Lab011	3	15.2	14.3	82	1.1	7.7	21.3	20	1.8		5.4	4.8	82	0.4	8.3	27.8	20	1.1	
Lab013	3	15.1	15.3	176	0.3	2.0	5.2	20	>6		5.4	5.5	176	0.2	3.6	9.1	20	5.0	
Lab014	3	14.4	14.3	244	0.3	2.1	4.9	20	>6		5.7	5.6	244	0.2	3.6	8.9	20	5.1	
Lab016	3	12.8	12.8	92	0.4	3.1	6.3	20	>6		4.2	4.3	92	0.2	4.7	11.7	20	3.7	
Lab029	3	14.8	15.1	123	0.6	4.0	10.0	20	4.5		5.1	5.3	123	0.3	5.7	15.2	20	2.8	
Lab044	3	14.7	14.7	83	0.5	3.4	6.8	20	5.9		5.1	5.1	83	0.2	3.9	7.8	20	5.1	
Lab051	3	14.4	15.4	84	2.0	13.0	32.9	20	1.0		5.7	5.3	84	1.5	28.3	63.6	20	0.5	
Total	-	-	14.6	1084	1.0	6.8	-	-	-		-	5.1	1084	0.6	11.8	-	-	-	

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%
總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)
σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

TOP

Reagent Kit	Reagent Code
Innovation	3

各檢驗單位 **G6PD** 定量檢驗 內部品管統計結果（總表二）

選擇 LotNo：BA1116N (2022-04-26 ~ 2100-12-31) Change

選擇 Reagent Kit：3 - Innovation Change

列印總表二

Control N 月統計與累計統計之比較表

		Control N (Lot No.: BA1116N)															
		Month (2025/06)								CUM (2022/06/01~2025/06/30)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab004	3	14.8	-	0	-	-	-	30	-	14.8	14.8	143	0.4	2.7	5.4	30	>6
Lab009	3	12.8	14.3	4	0.3	2.1	15.9	20	3.9	12.8	14.1	57	0.4	2.8	15.8	20	3.5
Lab011	3	15.2	-	0	-	-	-	20	-	15.2	14.3	82	1.1	7.7	21.3	20	1.8
Lab013	3	15.1	-	0	-	-	-	20	-	15.1	15.3	176	0.3	2.0	5.2	20	>6
Lab014	3	14.4	-	0	-	-	-	20	-	14.4	14.3	244	0.3	2.1	4.9	20	>6
Lab016	3	12.8	-	0	-	-	-	20	-	12.8	12.8	92	0.4	3.1	6.3	20	>6
Lab029	3	14.8	-	0	-	-	-	20	-	14.8	15.1	123	0.6	4.0	10.0	20	4.5
Lab044	3	14.7	15.2	4	0.1	0.7	4.7	20	>6	14.7	14.7	83	0.5	3.4	6.8	20	5.9
Lab051	3	14.4	17.8	3	0.3	1.7	27.0	20	-2.1	14.4	15.4	84	2.0	13.0	32.9	20	1.0
Total	-	-	15.6	11	1.5	9.6	-	-	-	-	14.6	1084	1.0	6.8	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計：2025 06 Change 累計期間 (CUM)：自 2022 06 01 至 2025 06 30 Change

[\[TOP\]](#)

Control D 月統計與累計統計之比較表

		Control D (Lot No.: BA1116D)															
		Month (2025/06)								CUM (2022/06/01~2025/06/30)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab004	3	4.6	-	0	-	-	-	30	-	4.6	4.6	143	0.2	4.3	8.7	30	>6
Lab009	3	4.2	5.0	4	0.1	2.0	23.0	20	0.5	4.2	5.1	57	0.2	3.9	29.3	20	-0.4
Lab011	3	5.4	-	0	-	-	-	20	-	5.4	4.8	82	0.4	8.3	27.8	20	1.1
Lab013	3	5.4	-	0	-	-	-	20	-	5.4	5.5	176	0.2	3.6	9.1	20	5.0
Lab014	3	5.7	-	0	-	-	-	20	-	5.7	5.6	244	0.2	3.6	8.9	20	5.1
Lab016	3	4.2	-	0	-	-	-	20	-	4.2	4.3	92	0.2	4.7	11.7	20	3.7
Lab029	3	5.1	-	0	-	-	-	20	-	5.1	5.3	123	0.3	5.7	15.2	20	2.8
Lab044	3	5.1	5.2	4	0.1	1.9	5.8	20	>6	5.1	5.1	83	0.2	3.9	7.8	20	5.1
Lab051	3	5.7	4.8	3	0.2	4.2	24.1	20	1.0	5.7	5.3	84	1.5	28.3	63.6	20	0.5
Total	-	-	5.0	11	0.2	4.0	-	-	-	-	5.1	1084	0.6	11.8	-	-	-

誤差 (Bias)% = { [| 累計平均值 (Mean) - 目標值 (Target) |] / 目標值 (Target) } x 100%

總誤差 (TE, Total Error)% = 誤差 (Bias)% + 2 x 累計變異係數 (CV)

σ (Sigma) = (TEa% - 誤差 (Bias)%) / 累計變異係數 (CV)

月統計：2025 06 Change 累計期間 (CUM)：自 2022 06 01 至 2025 06 30 Change

[\[TOP\]](#)

Reagent Kit	Reagent Code
Innovation	3