

QC Chart of Internal Quality Control (IQC)for G6PD Quantitative Test

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Print Table

Lab049

	Control N		Control D	
QC Control Lot No.	BJ0922N		BJ0922D	
Duration of the Analyzing	Month (2025/10)	CUM (2025/04/30~2025/10/30)	Month (2025/10)	CUM (2025/04/30~2025/10/30)
Runs (N)	4	26	4	26
Mean (U/gHb)	13.4	14.4	4.8	4.5
SD	0.3	0.6	0.3	0.2
CV (%)	2.2	4.2	6.3	4.4
Target Value (U/gHb)	13.8	13.8	4.7	4.7
Total Error (%)	7.4	12.7	14.6	13.1
TEa (%)	20	20	20	20
σ	>6	3.7	2.8	3.6

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2025 10 Change ; Cumulative : from 2025 04 30 to 2025 10 30 Change

TOP

Control N and Control D SDI
QC Chart

Lot No.: BJ0922N ; Duration : 2025-04-30 to 2025-10-30 ; Target : 13.8 ; SD : 1.38 (●)

Lot No.: BJ0922D ; Duration : 2025-04-30 to 2025-10-30 ; Target : 4.7 ; SD : 0.47 (□)

Lab049

Lab049

+3 S.D.

+2 S.D.

+1 S.D.

Target

-1 S.D.

-2 S.D.

-3 S.D.

Date

04-30 05-07 05-14 05-21 05-28 06-04 06-11 06-17 06-25 07-04 07-09 07-16 07-23 07-30 08-06 08-13 08-20 08-27 09-03 09-10 09-17 09-30 10-08 10-15 10-22 10-29

Month : 2025 10 Change ; Cumulative : from 2025 04 30 to 2025 10 30 Change

TOP

Control N SDI QC Chart

Lot No.: BJ0922N ; Duration : 2025-04-30 to 2025-10-30 ; Target : 13.8 ; SD : 1.38

Lab049

Lab049

+3 S.D.

+2 S.D.

+1 S.D.

Target

-1 S.D.

-2 S.D.

-3 S.D.

Date

04-30 05-07 05-14 05-21 05-28 06-04 06-11 06-17 06-25 07-04 07-09 07-16 07-23 07-30 08-06 08-13 08-20 08-27 09-03 09-10 09-17 09-30 10-08 10-15 10-22 10-29


Month : 2025 10 Change ; Cumulative : from 2025 04 30 to 2025 10 30 Change

TOP

Control D SDI QC Chart

Lot No.: BJ0922D ; Duration : 2025-04-30 to 2025-10-30 ; Target : 4.7 ; SD : 0.47

Lab049

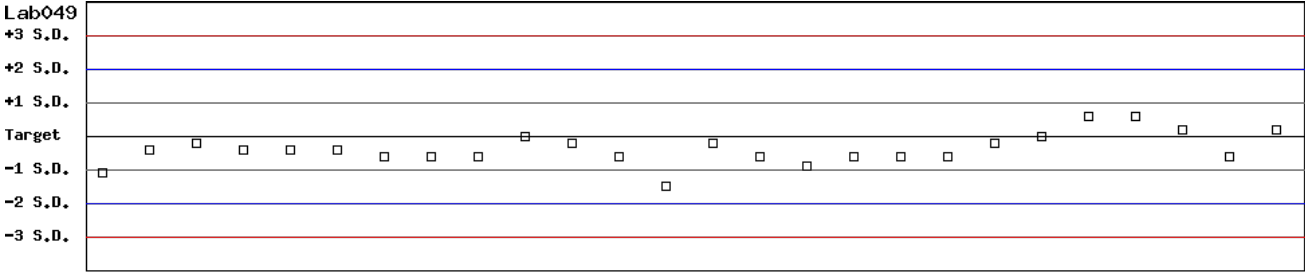

Date 04-30 05-07 05-14 05-21 05-28 06-04 06-11 06-17 06-25 07-04 07-09 07-16 07-23 07-30 08-06 08-13 08-20 08-27 09-03 09-10 09-17 09-30 10-08 10-15 10-22 10-29

Month : 2025 10 Change ; Cumulative : from 2025 04 30 to 2025 10 30 Change

[\[TOP\]](#)

Peer Group Statistics (Table 1)

Select LotNo : BJ0922N (2021-01-01 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

Print Table 1

Monthly

Month : 2025 10 Change

		Control N (Lot No.: BJ0922N)									Control D (Lot No.: BJ0922D)							
UnitID ↑	Reagent Kit (Code) ↑	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab034	5	13.8	15.3	3	0.1	0.7	12.2	20	>6		4.7	4.5	3	0.1	2.2	8.7	20	>6
Lab043	5	13.2	14.4	5	0.6	4.2	17.4	20	2.6		4.2	4.8	5	0.1	2.1	18.5	20	2.7
Lab049	5	13.8	13.4	4	0.3	2.2	7.4	20	>6		4.7	4.8	4	0.3	6.3	14.6	20	2.8
Lab050	5	13.8	14.0	3	0.3	2.1	5.7	20	>6		4.7	5.3	3	0.2	3.8	20.3	20	1.9
Total	-	-	14.2	15	0.8	5.6	-	-	-		-	4.8	15	0.3	6.3	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%
TE : Total Error(%) = Bias (%) + 2 x CV (%)
σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Cumulative

Cumulative : from 2021 01 01 to 2025 10 31 Change

		Control N (Lot No.: BJ0922N)									Control D (Lot No.: BJ0922D)							
UnitID ↑	Reagent Kit (Code) ↑	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	16.2	15.3	113	0.4	2.6	10.8	20	5.6		5.5	5.4	113	0.1	1.9	5.5	20	>6
Lab019	5	14.7	14.6	177	0.5	3.4	7.5	20	5.7		5.0	4.8	177	0.2	4.2	12.3	20	3.8
Lab020	5	14.4	14.4	90	0.5	3.5	6.9	20	5.7		4.6	4.7	90	0.2	4.3	10.7	20	4.1
Lab022	5	14.4	14.4	27	0.5	3.5	6.9	20	5.7		4.6	4.7	27	0.3	6.4	14.9	20	2.8
Lab025	5	14.4	14.2	13	0.5	3.5	8.4	20	5.3		5.4	5.1	13	0.3	5.9	17.3	20	2.4
Lab028	5	15.2	15.0	66	0.9	6.0	13.3	20	3.1		4.9	5.1	66	0.3	5.9	15.8	20	2.7
Lab032	5	14.9	14.3	153	0.4	2.8	9.6	20	5.7		4.8	4.5	153	0.2	4.4	15.1	20	3.1
Lab033	5	13.8	13.6	81	0.7	5.1	11.7	20	3.6		4.7	4.6	81	0.3	6.5	15.2	20	2.7
Lab034	5	13.8	14.5	49	0.6	4.1	13.3	20	3.6		4.7	5.0	49	0.7	14.0	34.4	20	1.0
Lab037	5	13.8	14.5	59	0.5	3.4	12.0	20	4.4		4.7	4.7	59	0.2	4.3	8.5	20	4.7
Lab040	5	14.4	14.8	169	0.2	1.4	5.5	20	>6		4.6	4.6	169	0.1	2.2	4.3	20	>6
Lab043	5	13.2	14.1	214	1.0	7.1	21.0	20	1.9		4.2	4.7	214	0.4	8.5	28.9	20	1.0
Lab049	5	13.8	14.4	281	0.6	4.2	12.7	20	3.7		4.7	4.7	281	0.3	6.4	12.8	20	3.1
Lab050	5	13.8	14.4	116	0.6	4.2	12.7	20	3.7		4.7	5.1	116	0.5	9.8	28.1	20	1.2
Total	-	-	14.5	1608	0.7	4.8	-	-	-		-	4.8	1608	0.4	8.3	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%
TE : Total Error(%) = Bias (%) + 2 x CV (%)
σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Medicon	5

Peer Group Statistics (Table 2)

Select LotNo :
BJ0922N (2021-01-01 ~ 2100-12-31)
Change

Select Reagent Kit :
5 - Medicon
Change

Print Table 2

Control N Month vs. Cumulative

		Control N (Lot No.: BJ0922N)															
		Month (2025/10)								CUM (2021/01/01~2025/10/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	16.2	-	0	-	-	-	20	-	16.2	15.3	113	0.4	2.6	10.8	20	5.6
Lab019	5	14.7	-	0	-	-	-	20	-	14.7	14.6	177	0.5	3.4	7.5	20	5.7
Lab020	5	14.4	-	0	-	-	-	20	-	14.4	14.4	90	0.5	3.5	6.9	20	5.7
Lab022	5	14.4	-	0	-	-	-	20	-	14.4	14.4	27	0.5	3.5	6.9	20	5.7
Lab025	5	14.4	-	0	-	-	-	20	-	14.4	14.2	13	0.5	3.5	8.4	20	5.3
Lab028	5	15.2	-	0	-	-	-	20	-	15.2	15.0	66	0.9	6.0	13.3	20	3.1
Lab032	5	14.9	-	0	-	-	-	20	-	14.9	14.3	153	0.4	2.8	9.6	20	5.7
Lab033	5	13.8	-	0	-	-	-	20	-	13.8	13.6	81	0.7	5.1	11.7	20	3.6
Lab034	5	13.8	15.3	3	0.1	0.7	12.2	20	>6	13.8	14.5	49	0.6	4.1	13.3	20	3.6
Lab037	5	13.8	-	0	-	-	-	20	-	13.8	14.5	59	0.5	3.4	12.0	20	4.4
Lab040	5	14.4	-	0	-	-	-	20	-	14.4	14.8	169	0.2	1.4	5.5	20	>6
Lab043	5	13.2	14.4	5	0.6	4.2	17.4	20	2.6	13.2	14.1	214	1.0	7.1	21.0	20	1.9
Lab049	5	13.8	13.4	4	0.3	2.2	7.4	20	>6	13.8	14.4	281	0.6	4.2	12.7	20	3.7
Lab050	5	13.8	14.0	3	0.3	2.1	5.7	20	>6	13.8	14.4	116	0.6	4.2	12.7	20	3.7
Total	-	-	14.2	15	0.8	5.6	-	-	-	-	14.5	1608	0.7	4.8	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month :
2025
10
Change

Cumulative : from
2021
01
01
to
2025
10
31
Change

[TOP](#)

Control D Month vs. Cumulative

		Control D (Lot No.: BJ0922D)															
		Month (2025/10)								CUM (2021/01/01~2025/10/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	5.5	-	0	-	-	-	20	-	5.5	5.4	113	0.1	1.9	5.5	20	>6
Lab019	5	5.0	-	0	-	-	-	20	-	5.0	4.8	177	0.2	4.2	12.3	20	3.8
Lab020	5	4.6	-	0	-	-	-	20	-	4.6	4.7	90	0.2	4.3	10.7	20	4.1
Lab022	5	4.6	-	0	-	-	-	20	-	4.6	4.7	27	0.3	6.4	14.9	20	2.8
Lab025	5	5.4	-	0	-	-	-	20	-	5.4	5.1	13	0.3	5.9	17.3	20	2.4
Lab028	5	4.9	-	0	-	-	-	20	-	4.9	5.1	66	0.3	5.9	15.8	20	2.7
Lab032	5	4.8	-	0	-	-	-	20	-	4.8	4.5	153	0.2	4.4	15.1	20	3.1
Lab033	5	4.7	-	0	-	-	-	20	-	4.7	4.6	81	0.3	6.5	15.2	20	2.7
Lab034	5	4.7	4.5	3	0.1	2.2	8.7	20	>6	4.7	5.0	49	0.7	14.0	34.4	20	1.0
Lab037	5	4.7	-	0	-	-	-	20	-	4.7	4.7	59	0.2	4.3	8.5	20	4.7
Lab040	5	4.6	-	0	-	-	-	20	-	4.6	4.6	169	0.1	2.2	4.3	20	>6
Lab043	5	4.2	4.8	5	0.1	2.1	18.5	20	2.7	4.2	4.7	214	0.4	8.5	28.9	20	1.0
Lab049	5	4.7	4.8	4	0.3	6.3	14.6	20	2.8	4.7	4.7	281	0.3	6.4	12.8	20	3.1
Lab050	5	4.7	5.3	3	0.2	3.8	20.3	20	1.9	4.7	5.1	116	0.5	9.8	28.1	20	1.2
Total	-	-	4.8	15	0.3	6.3	-	-	-	-	4.8	1608	0.4	8.3	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month :
2025
10
Change

Cumulative : from
2021
01
01
to
2025
10
31
Change

[TOP](#)

Reagent Kit	Reagent Code
Medicon	5