

QC Chart of Internal Quality Control (IQC)for G6PD Quantitative Test

Select LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

Print Table

Lab037

QC Control Lot No.	Control N		Control D	
	BA1116N		BA1116D	
	Month (2025/12)	CUM (2024/11/08~2025/12/30)	Month (2025/12)	CUM (2024/11/08~2025/12/30)
Duration of the Analyzing				
Runs (N)	2	26	2	26
Mean (U/gHb)	13.9	13.7	4.3	4.4
SD	-	0.3	-	0.2
CV (%)	-	2.2	-	4.5
Target Value (U/gHb)	12.8	12.8	4.2	4.2
Total Error (%)	-	11.4	-	13.9
TEa (%)	20	20	20	20
σ	-	5.9	-	3.4

Bias (%) = [(| Mean - Target |) / Target] x 100%
TE : Total Error(%) = Bias (%) + 2 × CV (%)
σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2025 12 Change ; Cumulative : from 2024 11 08 to 2025 12 30 Change


TOP

Control N and Control D SDI
QC Chart

Lot No.: BA1116N ; Duration : 2024-11-08 to 2025-12-30 ; Target : 12.8 ; SD : 1.28 (●)

Lot No.: BA1116D ; Duration : 2024-11-08 to 2025-12-30 ; Target : 4.2 ; SD : 0.42 (□)

Lab037

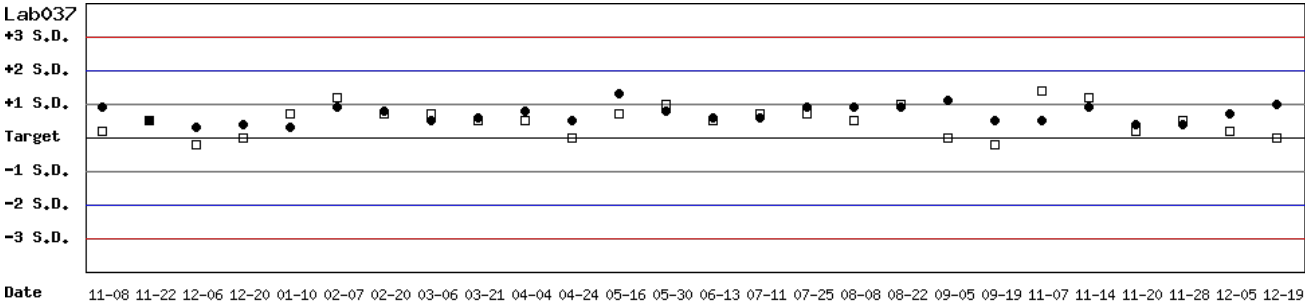


Month : 2025 12 Change ; Cumulative : from 2024 11 08 to 2025 12 30 Change

TOP

Control N SDI QC Chart

Lot No.: BA1116N ; Duration : 2024-11-08 to 2025-12-30 ; Target : 12.8 ; SD : 1.28

Lab037

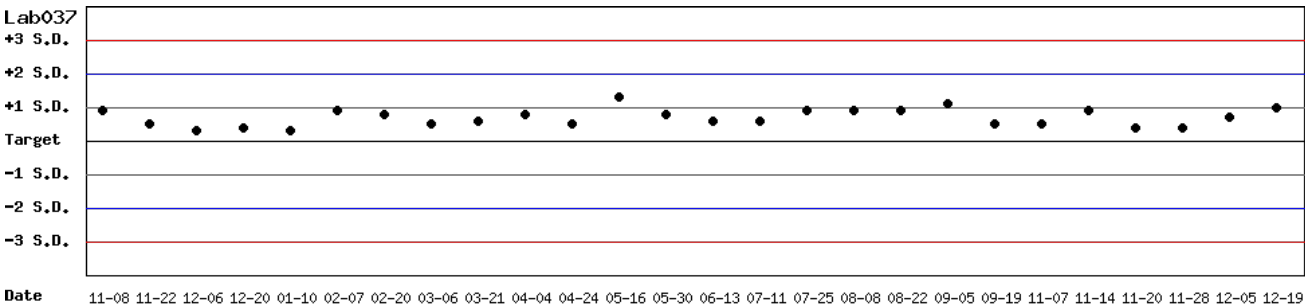


Month : 2025 12 Change ; Cumulative : from 2024 11 08 to 2025 12 30 Change

TOP

Control D SDI QC Chart

Lot No.: BA1116D ; Duration : 2024-11-08 to 2025-12-30 ; Target : 4.2 ; SD : 0.42

Lab037

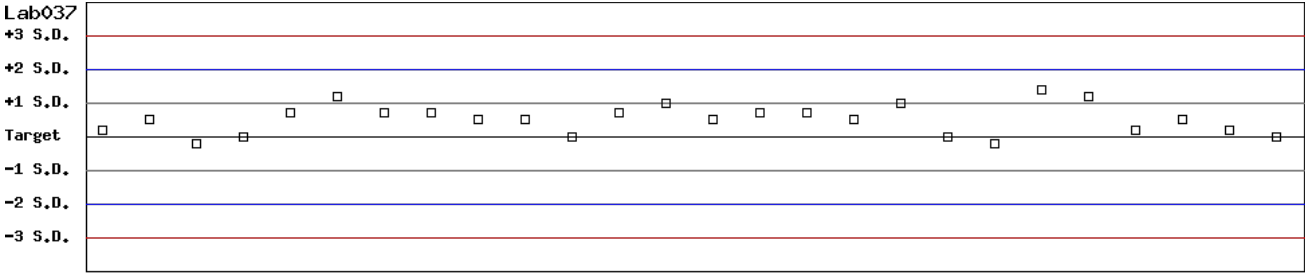

Date 11-08 11-22 12-06 12-20 01-10 02-07 02-20 03-06 03-21 04-04 04-24 05-16 05-30 06-13 07-11 07-25 08-08 08-22 09-05 09-19 11-07 11-14 11-20 11-28 12-05 12-19

Month : 2025 12 Change ; Cumulative : from 2024 11 08 to 2025 12 30 Change

[TOP](#)

Peer Group Statistics (Table 1)

Select LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

[Print Table 1](#)

Monthly

Month : 2025 12 Change

		Control N (Lot No.: BA1116N)								Control D (Lot No.: BA1116D)							
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab037	5	12.8	13.9	2	-	-	-	20	-	4.2	4.3	2	-	-	-	20	-
Total	-	-	13.9	2	-	-	-	-	-	-	4.3	2	-	-	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Cumulative

Cumulative : from 2022 06 01 to 2025 12 31 Change

		Control N (Lot No.: BA1116N)								Control D (Lot No.: BA1116D)							
UnitID	Reagent Kit (Code)	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	13.6	13.9	147	0.3	2.2	6.5	20	>6	4.3	4.4	147	0.1	2.3	6.9	20	>6
Lab027	5	12.8	14.9	74	0.9	6.0	28.5	20	0.6	4.2	4.9	74	0.4	8.2	33.0	20	0.4
Lab032	5	14.9	13.8	39	0.4	2.9	13.2	20	4.4	5.1	4.3	39	0.2	4.7	25.0	20	0.9
Lab033	5	15.8	15.7	85	0.6	3.8	8.3	20	5.1	4.7	4.8	85	0.2	4.2	10.5	20	4.3
Lab037	5	12.8	13.8	54	0.4	2.9	13.6	20	4.2	4.2	4.4	54	0.2	4.5	13.9	20	3.4
Lab052	5	15.3	15.3	15	0.3	2.0	3.9	20	>6	4.6	4.7	15	0.2	4.3	10.7	20	4.1
Total	-	-	14.5	418	0.9	6.2	-	-	-	-	4.6	418	0.3	6.5	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

[\[TOP\]](#)

Reagent Kit	Reagent Code
Medicon	5

Peer Group Statistics (Table 2)

Select LotNo : BA1116N (2022-04-26 ~ 2100-12-31) Change

Select Reagent Kit : 5 - Medicon Change

Print Table 2

Control N Month vs. Cumulative

		Control N (Lot No.: BA1116N)																
		Month (2025/12)									CUM (2022/06/01~2025/12/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	13.6	-	0	-	-	-	20	-		13.6	13.9	147	0.3	2.2	6.5	20	>6
Lab027	5	12.8	-	0	-	-	-	20	-		12.8	14.9	74	0.9	6.0	28.5	20	0.6
Lab032	5	14.9	-	0	-	-	-	20	-		14.9	13.8	39	0.4	2.9	13.2	20	4.4
Lab033	5	15.8	-	0	-	-	-	20	-		15.8	15.7	85	0.6	3.8	8.3	20	5.1
Lab037	5	12.8	13.9	2	-	-	-	20	-		12.8	13.8	54	0.4	2.9	13.6	20	4.2
Lab052	5	15.3	-	0	-	-	-	20	-		15.3	15.3	15	0.3	2.0	3.9	20	>6
Total	-	-	13.9	2	-	-	-	-	-		-	14.5	418	0.9	6.2	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2025 12 Change Cumulative : from 2022 06 01 to 2025 12 31 Change

TOP

Control D Month vs. Cumulative

		Control D (Lot No.: BA1116D)																
		Month (2025/12)									CUM (2022/06/01~2025/12/31)							
UnitID 	Reagent Kit (Code) 	Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ		Target (U/gHb)	Mean (U/gHb)	n for Mean	SD	CV (%)	TE (%)	TEa (%)	σ
Lab018	5	4.3	-	0	-	-	-	20	-		4.3	4.4	147	0.1	2.3	6.9	20	>6
Lab027	5	4.2	-	0	-	-	-	20	-		4.2	4.9	74	0.4	8.2	33.0	20	0.4
Lab032	5	5.1	-	0	-	-	-	20	-		5.1	4.3	39	0.2	4.7	25.0	20	0.9
Lab033	5	4.7	-	0	-	-	-	20	-		4.7	4.8	85	0.2	4.2	10.5	20	4.3
Lab037	5	4.2	4.3	2	-	-	-	20	-		4.2	4.4	54	0.2	4.5	13.9	20	3.4
Lab052	5	4.6	-	0	-	-	-	20	-		4.6	4.7	15	0.2	4.3	10.7	20	4.1
Total	-	-	4.3	2	-	-	-	-	-		-	4.6	418	0.3	6.5	-	-	-

Bias (%) = [(| Mean - Target |) / Target] x 100%

TE : Total Error(%) = Bias (%) + 2 × CV (%)

σ (Sigma) = [TEa% - Bias (%)] / CV (%)

Month : 2025 12 Change Cumulative : from 2022 06 01 to 2025 12 31 Change

TOP

Reagent Kit	Reagent Code
Medicon	5